

Metadata: Winston-Lutz Analysis

Função: Físico Number of images: 8

Acelerador Linear: Synergy Full Maximum 2D CAX->BB distance: 1.91mm Modelo MLC: Agility Median 2D CAX->BB distance: 1.03mm

Shift to iso: facing gantry, move BB: RIGHT 0.29mm; IN 0.08mm; UP 0.47mm

Gantry 3D isocenter diameter: 1.20mm (4/8 images considered)

Maximum Gantry RMS deviation (mm): 0.87mm Maximum EPID RMS deviation (mm): 2.53mm


Gantry+Collimator 3D isocenter diameter: 1.37mm (6/8 images considered)

Collimator 2D isocenter diameter: 0.47mm (3/8 images considered)

Maximum Collimator RMS deviation (mm): 1.09

Couch 2D isocenter diameter: 1.61mm (3/8 images considered)

Maximum Couch RMS deviation (mm): 1.91

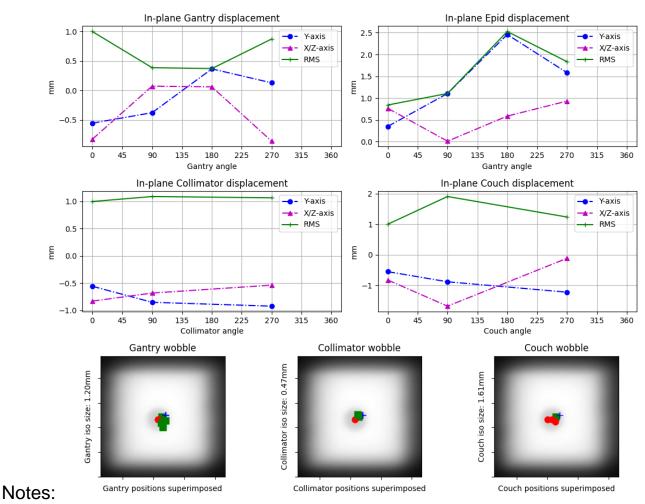

Análise de Winston-Lutz com pylinac

Campo 2 cm X 2 cm Tolerância: 1 mm Limite de ação: 0.8 mm

Local das imagens: ./WL\WL_2020-06-01_analisado




Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

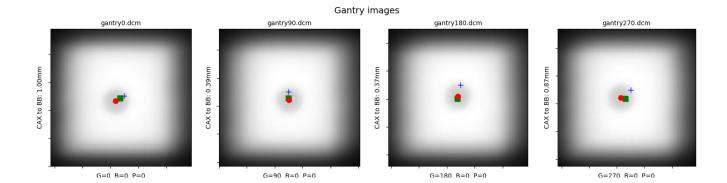




Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

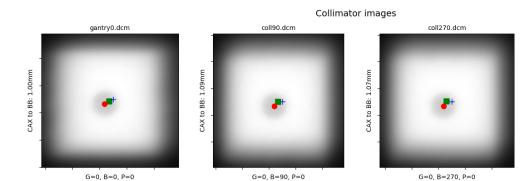



Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

