



# Winston-Lutz Analysis

## Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

## Winston-Lutz Analysis

Number of images: 8

Maximum 2D CAX->BB distance: 1.42mm

Median 2D CAX->BB distance: 0.75mm

Shift to iso: facing gantry, move BB: RIGHT 0.47mm; OUT 0.15mm; DOWN 0.47mm

Gantry 3D isocenter diameter: 1.25mm (4/8 images considered)

Maximum Gantry RMS deviation (mm): 1.14mm

Maximum EPID RMS deviation (mm): 3.03mm






Gantry+Collimator 3D isocenter diameter: 1.82mm (6/8 images considered)

Collimator 2D isocenter diameter: 0.57mm (3/8 images considered)


Maximum Collimator RMS deviation (mm): 0.96

Couch 2D isocenter diameter: 1.52mm (3/8 images considered)


Maximum Couch RMS deviation (mm): 1.42

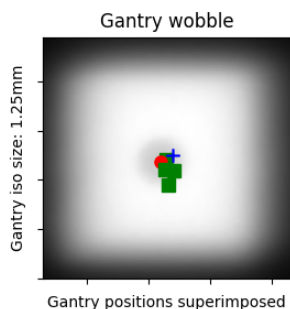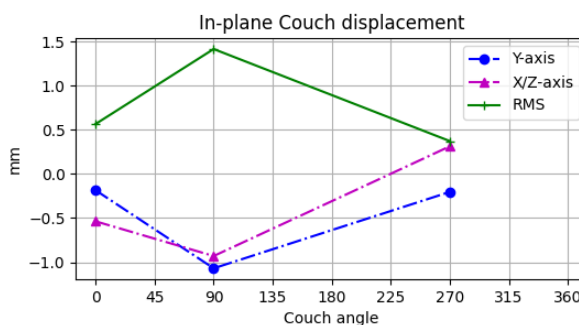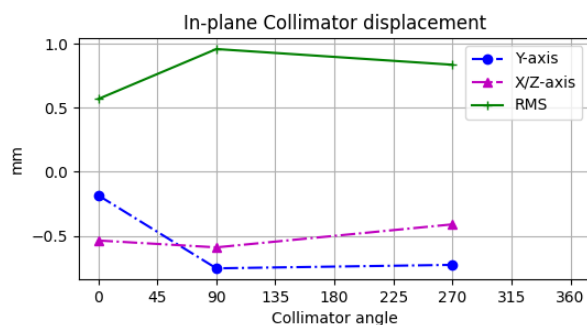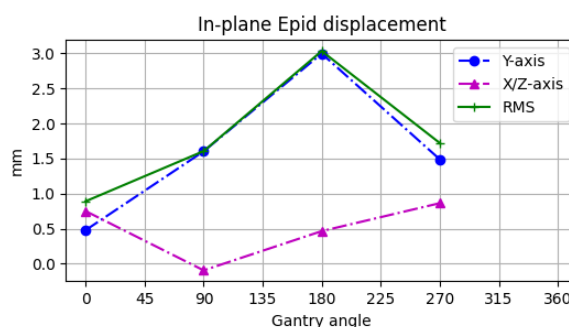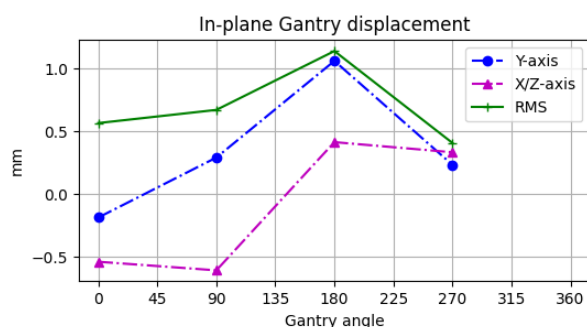

Gantry positions superimposed

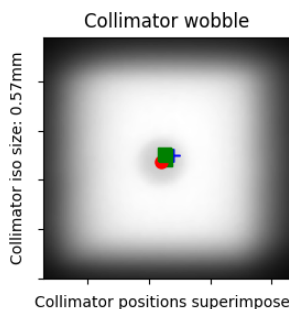

Collimator positions superimposed

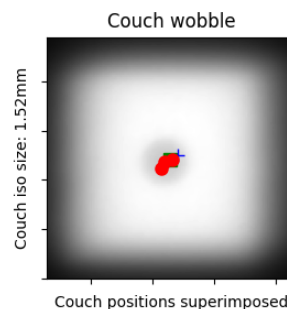

Couch positions superimposed

## Notes:

Análise de Winston-Lutz com pylinac

Campo 2 cm X 2 cm

Tolerância: 1 mm

Limite de ação: 0.8 mm

Local das imagens: ./WL\WL\_2020-08-17



# Winston-Lutz Analysis

## Metadata:


Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

Gantry images

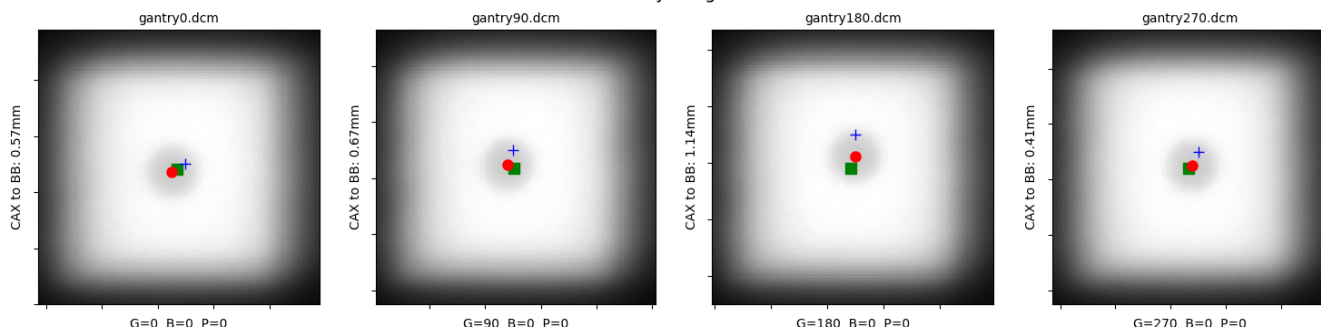



# Winston-Lutz Analysis

## Metadata:


Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

Collimator images

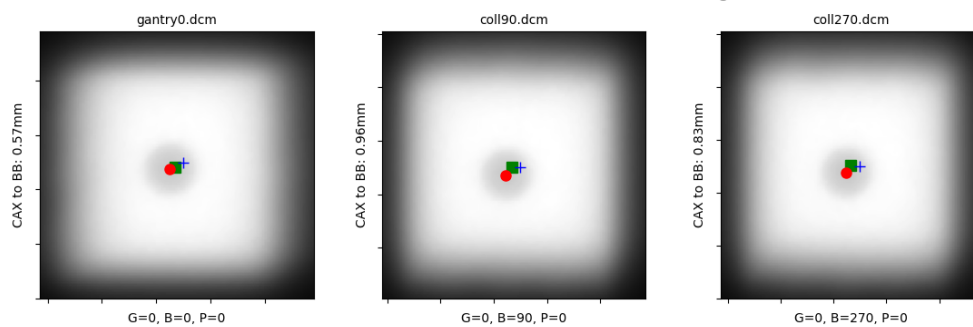



# Winston-Lutz Analysis

## Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

Couch images

