

Metadata: Winston-Lutz Analysis

Função: Físico Number of images: 8

Acelerador Linear: Synergy Full Maximum 2D CAX->BB distance: 1.74mm Modelo MLC: Agility Median 2D CAX->BB distance: 0.83mm

Shift to iso: facing gantry, move BB: LEFT 0.13mm; OUT 0.10mm; DOWN 0.27mm

Gantry 3D isocenter diameter: 1.17mm (4/8 images considered)

Maximum Gantry RMS deviation (mm): 0.91mm Maximum EPID RMS deviation (mm): 2.91mm


Gantry+Collimator 3D isocenter diameter: 1.75mm (6/8 images considered)

Collimator 2D isocenter diameter: 0.58mm (3/8 images considered)

Maximum Collimator RMS deviation (mm): 0.84

Couch 2D isocenter diameter: 1.79mm (3/8 images considered)

Maximum Couch RMS deviation (mm): 1.74

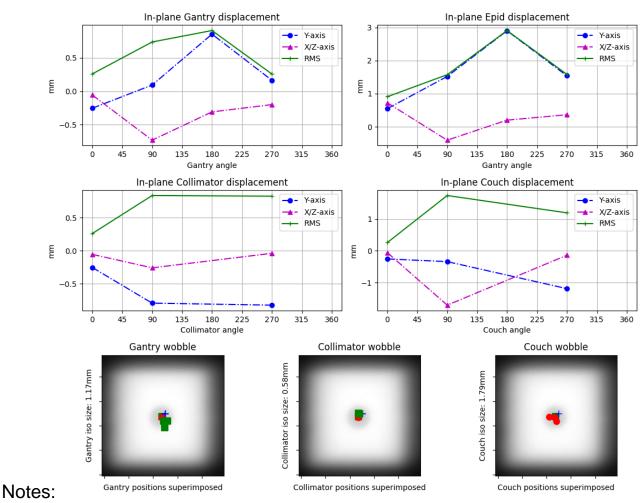

Análise de Winston-Lutz com pylinac

Campo 2 cm X 2 cm
Tolerância: 1 mm

Limite de ação: 0.8 mm

Local das imagens: ./WL\WL\_2020-05-04\_analisado




Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

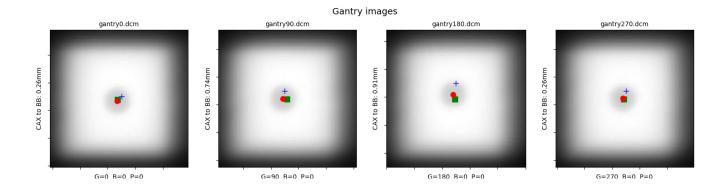




Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

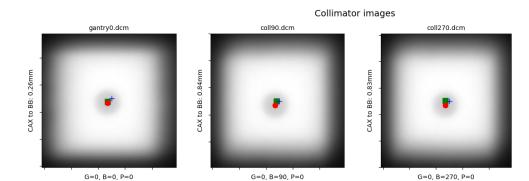



Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

