



# Winston-Lutz Analysis

## Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

## Winston-Lutz Analysis

Number of images: 8

Maximum 2D CAX->BB distance: 0.77mm

Median 2D CAX->BB distance: 0.60mm

Shift to iso: facing gantry, move BB: RIGHT 0.19mm; IN 0.34mm; DOWN 0.01mm

Gantry 3D isocenter diameter: 0.62mm (4/8 images considered)

Maximum Gantry RMS deviation (mm): 0.64mm

Maximum EPID RMS deviation (mm): 1.54mm






Gantry+Collimator 3D isocenter diameter: 0.80mm (6/8 images considered)

Collimator 2D isocenter diameter: 0.80mm (3/8 images considered)


Maximum Collimator RMS deviation (mm): 0.62

Couch 2D isocenter diameter: 0.69mm (3/8 images considered)


Maximum Couch RMS deviation (mm): 0.77

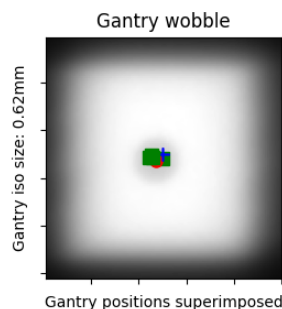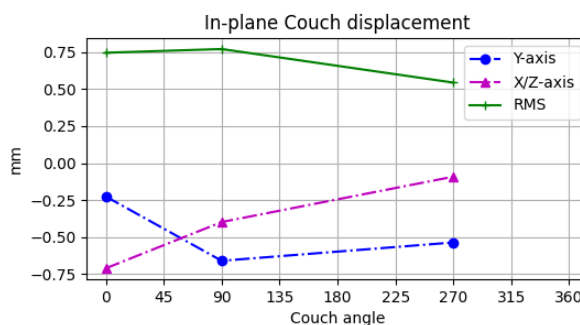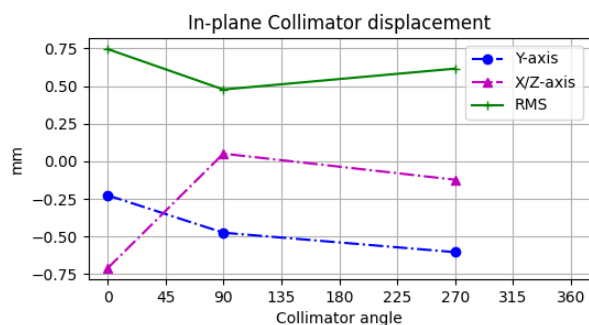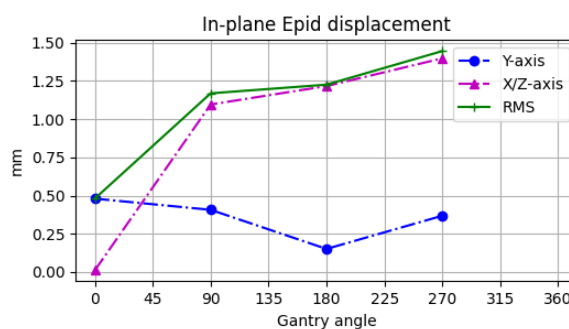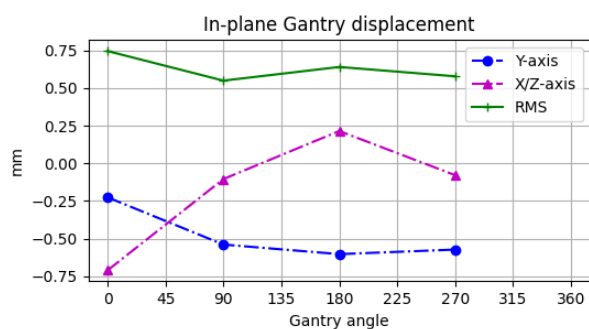

Gantry positions superimposed

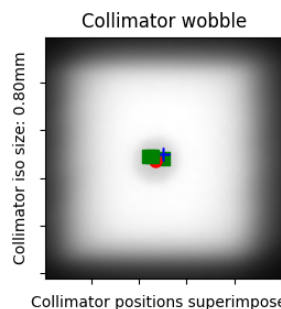

Collimator positions superimposed

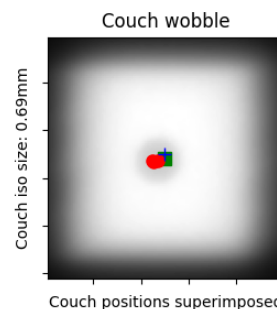

Couch positions superimposed

## Notes:

Análise de Winston-Lutz com pylinac

Campo 2 cm X 2 cm

Tolerância: 1 mm

Limite de ação: 0.8 mm

Local das imagens: ./WL\WL\_2020-07-24\_r



# Winston-Lutz Analysis

## Metadata:


Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

Gantry images

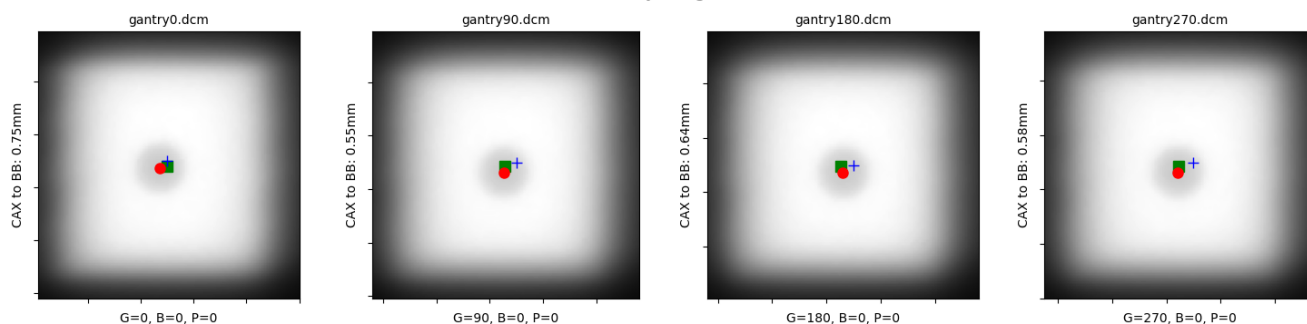



# Winston-Lutz Analysis

## Metadata:


Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

Collimator images

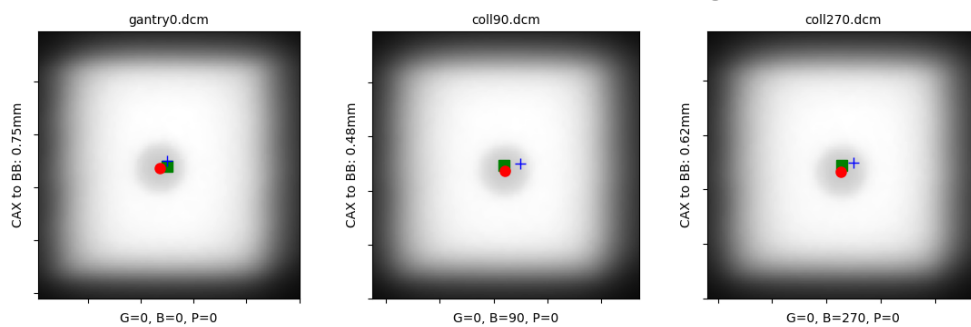



# Winston-Lutz Analysis

## Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

Couch images

