



Winston-Lutz Analysis

Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

Winston-Lutz Analysis

Number of images: 8

Maximum 2D CAX->BB distance: 1.42mm

Median 2D CAX->BB distance: 0.98mm

Shift to iso: facing gantry, move BB: RIGHT 0.06mm; IN 0.33mm; DOWN 0.39mm

Gantry 3D isocenter diameter: 1.14mm (4/8 images considered)

Maximum Gantry RMS deviation (mm): 0.68mm

Maximum EPID RMS deviation (mm): 7.42mm






Gantry+Collimator 3D isocenter diameter: 1.59mm (6/8 images considered)

Collimator 2D isocenter diameter: 0.45mm (3/8 images considered)


Maximum Collimator RMS deviation (mm): 1.27

Couch 2D isocenter diameter: 0.96mm (3/8 images considered)


Maximum Couch RMS deviation (mm): 1.42

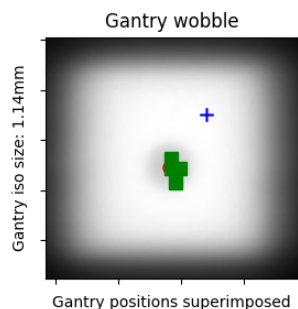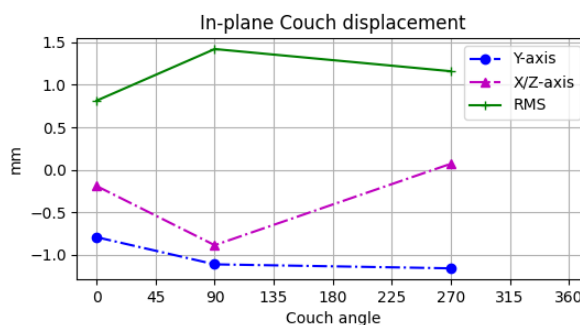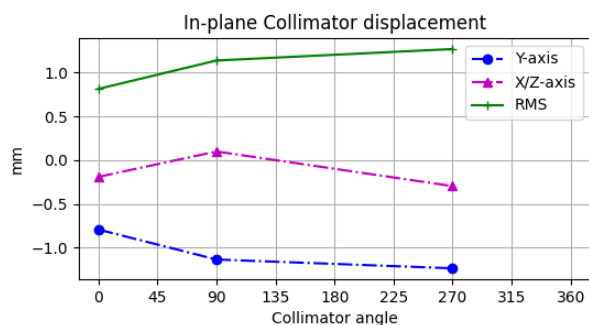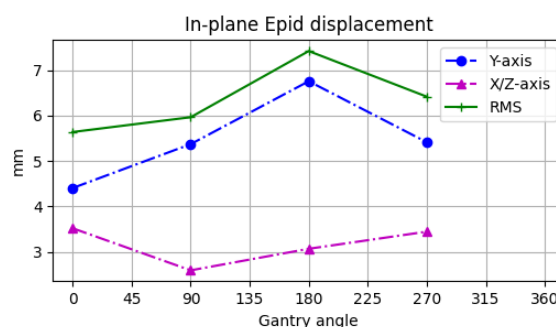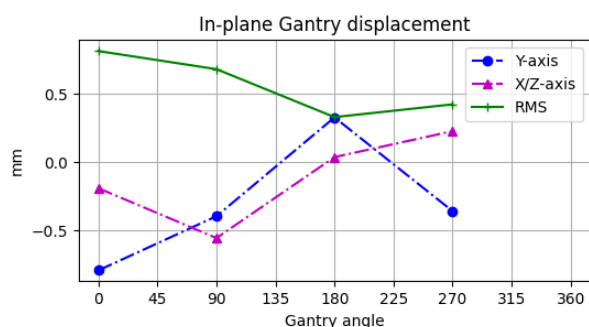

Gantry positions superimposed

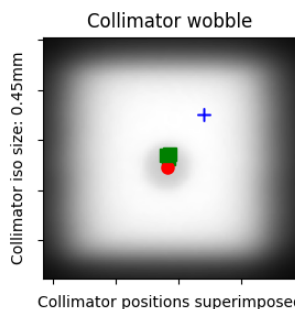

Collimator positions superimposed

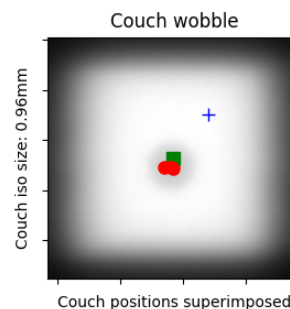

Couch positions superimposed

Notes:

Análise de Winston-Lutz com pylinac

Campo 2 cm X 2 cm

Tolerância: 1 mm

Limite de ação: 0.8 mm

Local das imagens: ./WL\WL_2020-01-14_analisado



Winston-Lutz Analysis

Metadata:


Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

Gantry images

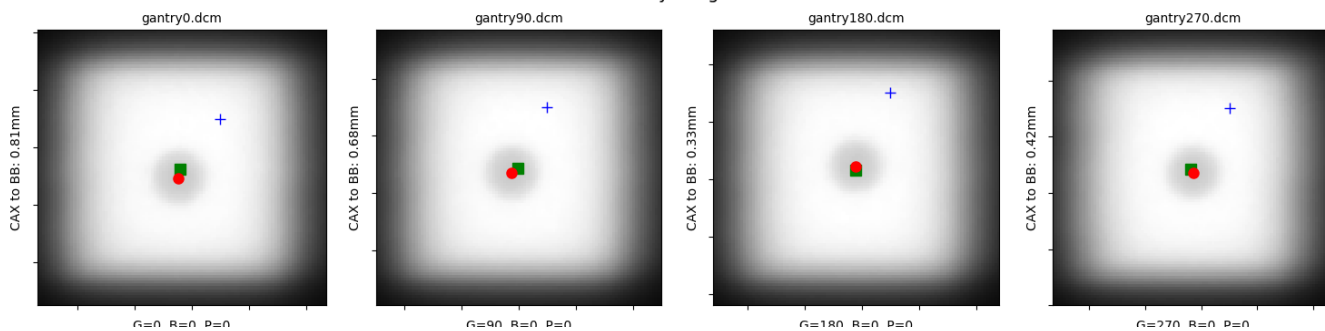



Winston-Lutz Analysis

Metadata:


Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

Collimator images

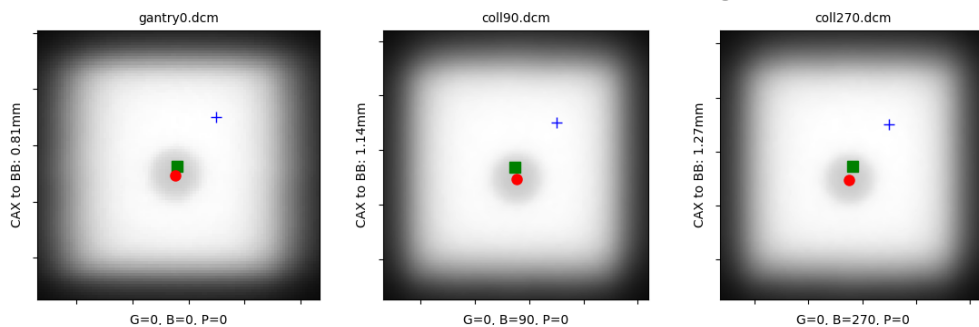



Winston-Lutz Analysis

Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

Couch images

