

Metadata: Winston-Lutz Analysis

Função: Físico Number of images: 9

Acelerador Linear: Synergy Full Maximum 2D CAX->BB distance: 1.11mm Modelo MLC: Agility Median 2D CAX->BB distance: 0.77mm

Shift to iso: facing gantry, move BB: LEFT 0.11mm; IN 0.24mm; DOWN 0.09mm

Gantry 3D isocenter diameter: 0.93mm (4/9 images considered)

Maximum Gantry RMS deviation (mm): 0.56mm Maximum EPID RMS deviation (mm): 2.54mm


Gantry+Collimator 3D isocenter diameter: 1.63mm (7/9 images considered)

Collimator 2D isocenter diameter: 0.73mm (4/9 images considered)

Maximum Collimator RMS deviation (mm): 1.11

Couch 2D isocenter diameter: 0.99mm (3/9 images considered)

Maximum Couch RMS deviation (mm): 0.80

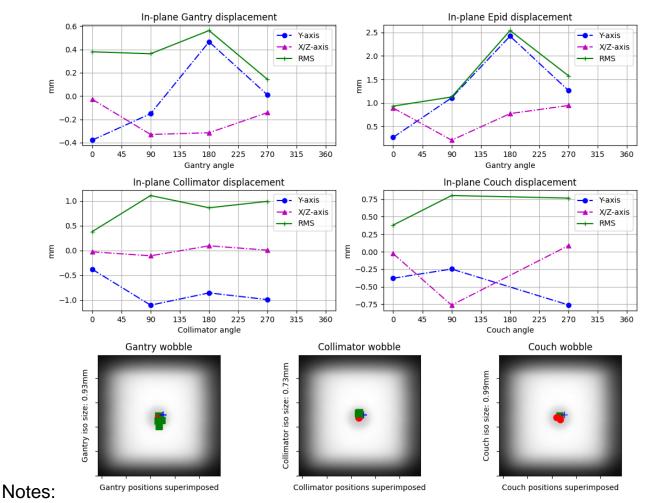

Análise de Winston-Lutz com pylinac

Campo 2 cm X 2 cm Tolerância: 1 mm Limite de ação: 0.8 mm

Local das imagens: ./WL\_IOSP\_20201216\WL\_2020-11-27\_tarde




Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

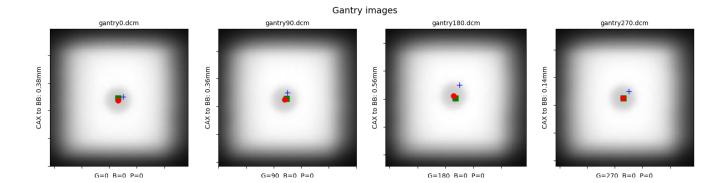




Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

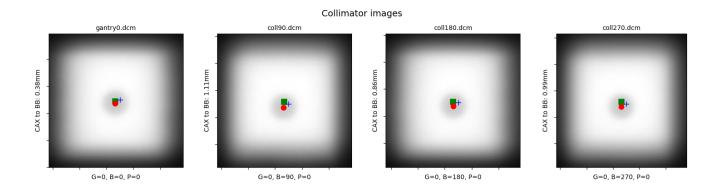



Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

