

Metadata: Winston-Lutz Analysis

Função: Físico Number of images: 9

Acelerador Linear: Synergy Full Maximum 2D CAX->BB distance: 1.47mm Modelo MLC: Agility Median 2D CAX->BB distance: 0.64mm

Shift to iso: facing gantry, move BB: RIGHT 0.27mm; OUT 0.15mm; DOWN 1.19mm

Gantry 3D isocenter diameter: 1.07mm (4/9 images considered)

Maximum Gantry RMS deviation (mm): 1.47mm Maximum EPID RMS deviation (mm): 3.10mm


Gantry+Collimator 3D isocenter diameter: 1.63mm (7/9 images considered)

Collimator 2D isocenter diameter: 0.56mm (4/9 images considered)

Maximum Collimator RMS deviation (mm): 0.64

Couch 2D isocenter diameter: 0.87mm (3/9 images considered)

Maximum Couch RMS deviation (mm): 0.91

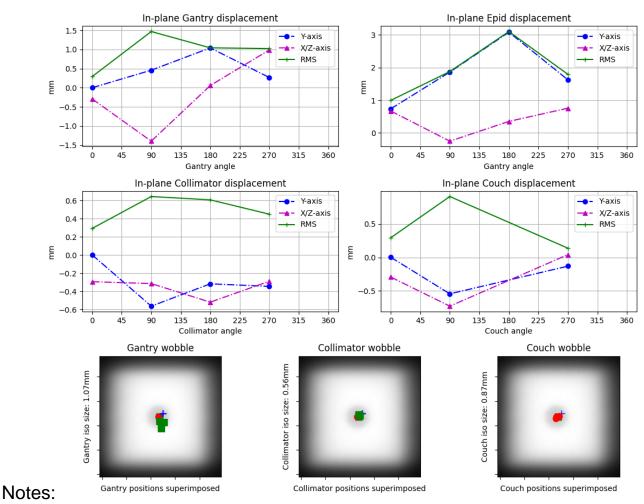

Análise de Winston-Lutz com pylinac

Campo 2 cm X 2 cm
Tolerância: 1 mm

Limite de ação: 0.8 mm

Local das imagens: ./WL_IOSP_20201216\WL_2020-10-29




Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

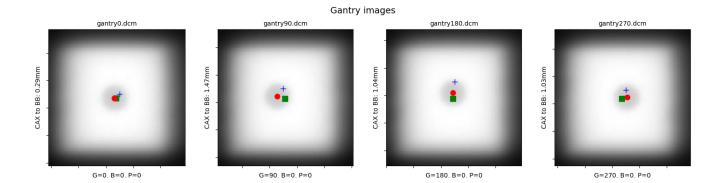




Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

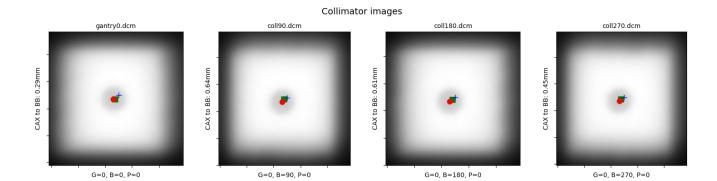



Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

