



Winston-Lutz Analysis

Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

Winston-Lutz Analysis

Number of images: 8

Maximum 2D CAX->BB distance: 1.26mm

Median 2D CAX->BB distance: 0.74mm

Shift to iso: facing gantry, move BB: RIGHT 0.10mm; OUT 0.05mm; UP 0.19mm

Gantry 3D isocenter diameter: 1.25mm (4/8 images considered)

Maximum Gantry RMS deviation (mm): 0.80mm

Maximum EPID RMS deviation (mm): 2.86mm








Gantry+Collimator 3D isocenter diameter: 1.58mm (6/8 images considered)

Collimator 2D isocenter diameter: 0.33mm (3/8 images considered)

Maximum Collimator RMS deviation (mm): 0.78

Couch 2D isocenter diameter: 1.46mm (3/8 images considered)

Maximum Couch RMS deviation (mm): 1.26

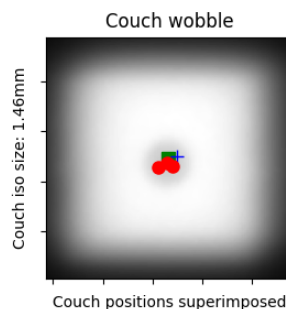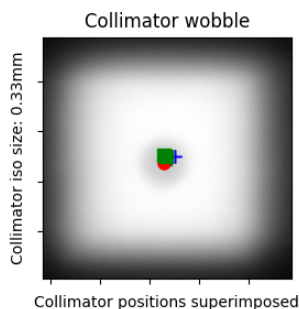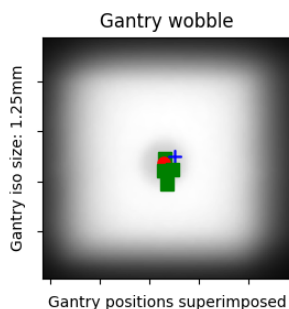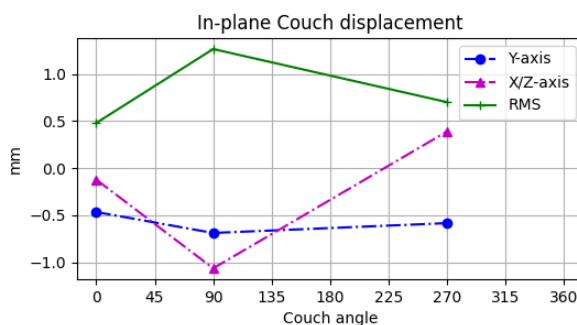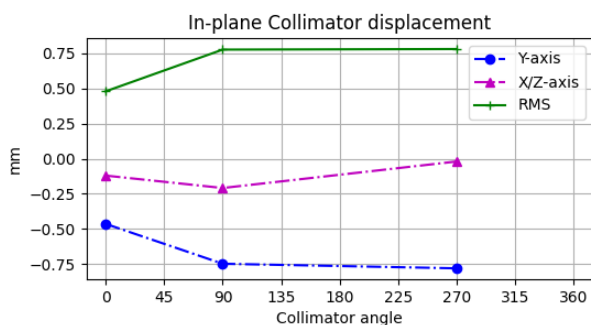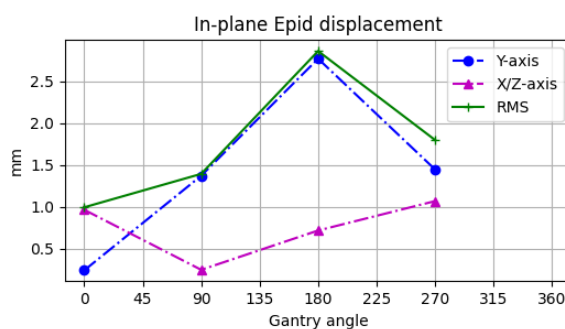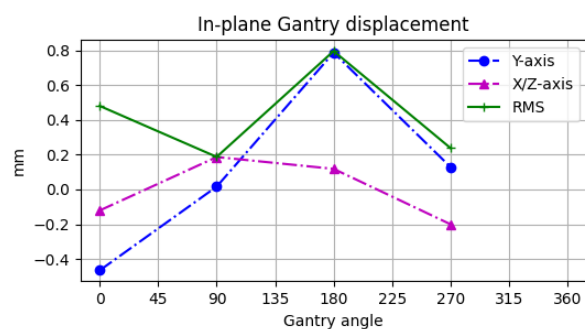

Notes:

Análise de Winston-Lutz com pylinac

Campo 2 cm X 2 cm

Tolerância: 1 mm

Limite de ação: 0.8 mm

Local das imagens: ./WL\WL_2020-08-14_analisado



Winston-Lutz Analysis

Metadata:


Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

Gantry images

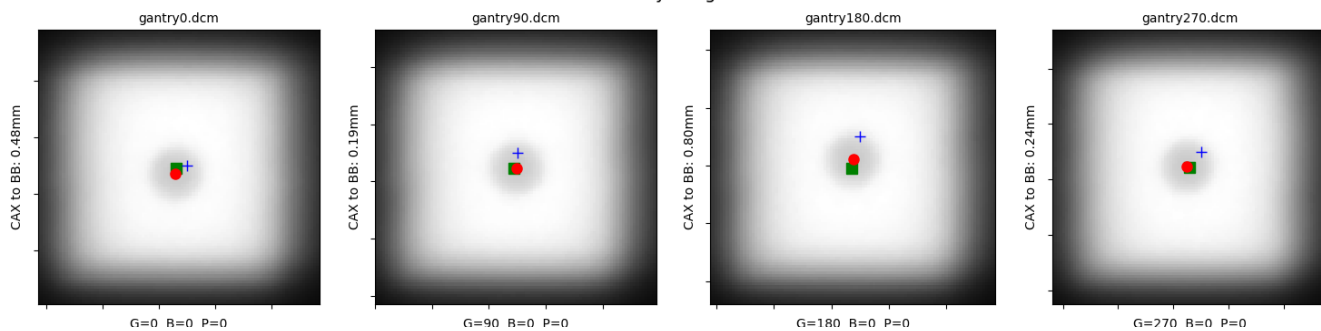



Winston-Lutz Analysis

Metadata:


Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

Collimator images

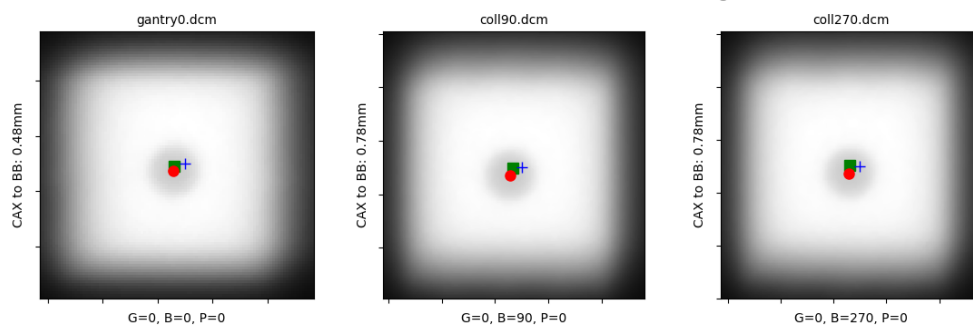



Winston-Lutz Analysis

Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

Couch images

