

Metadata: Winston-Lutz Analysis

Função: Físico Number of images: 9

Acelerador Linear: Synergy Full Maximum 2D CAX->BB distance: 0.68mm Modelo MLC: Agility Median 2D CAX->BB distance: 0.51mm

Shift to iso: facing gantry, move BB: RIGHT 0.08mm; IN 0.13mm; UP 0.20mm

Gantry 3D isocenter diameter: 1.11mm (4/9 images considered)

Maximum Gantry RMS deviation (mm): 0.51mm Maximum EPID RMS deviation (mm): 2.38mm


Gantry+Collimator 3D isocenter diameter: 1.17mm (7/9 images considered)

Collimator 2D isocenter diameter: 0.29mm (4/9 images considered)

Maximum Collimator RMS deviation (mm): 0.68

Couch 2D isocenter diameter: 0.60mm (3/9 images considered)

Maximum Couch RMS deviation (mm): 0.62

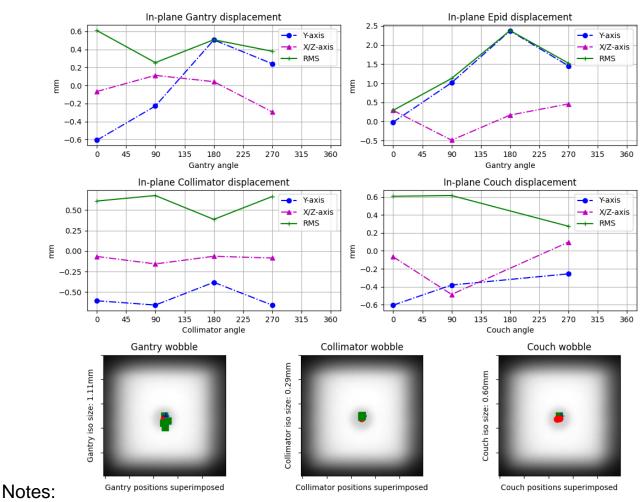

Análise de Winston-Lutz com pylinac

Campo 2 cm X 2 cm Tolerância: 1 mm Limite de ação: 0.8 mm

Local das imagens: ./WL3\WL\_2020-09-10\_r




Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

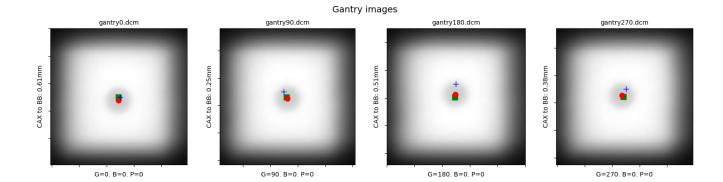




Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

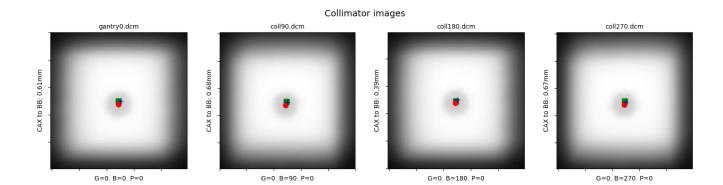



Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

