



Winston-Lutz Analysis

Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

Winston-Lutz Analysis

Number of images: 8

Maximum 2D CAX->BB distance: 1.31mm

Median 2D CAX->BB distance: 1.00mm

Shift to iso: facing gantry, move BB: LEFT 0.34mm; IN 0.27mm; DOWN 0.93mm

Gantry 3D isocenter diameter: 1.08mm (4/8 images considered)

Maximum Gantry RMS deviation (mm): 1.16mm

Maximum EPID RMS deviation (mm): 2.63mm








Gantry+Collimator 3D isocenter diameter: 1.56mm (6/8 images considered)

Collimator 2D isocenter diameter: 0.49mm (3/8 images considered)

Maximum Collimator RMS deviation (mm): 1.11

Couch 2D isocenter diameter: 1.47mm (3/8 images considered)

Maximum Couch RMS deviation (mm): 1.31

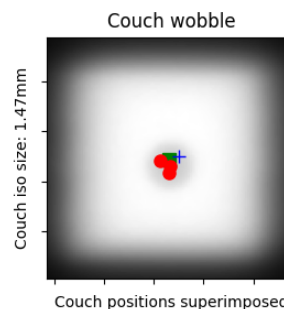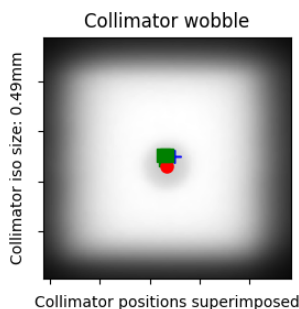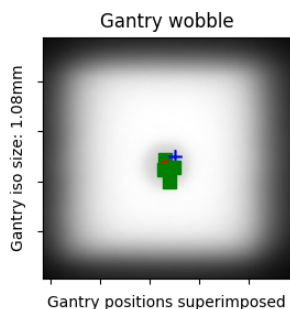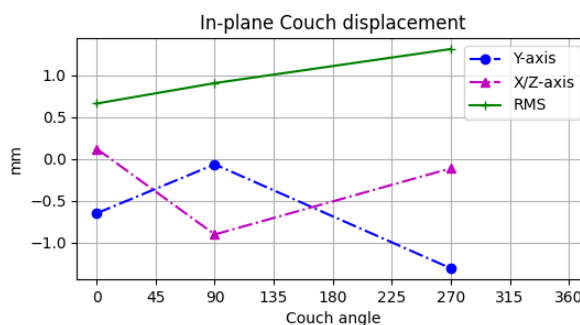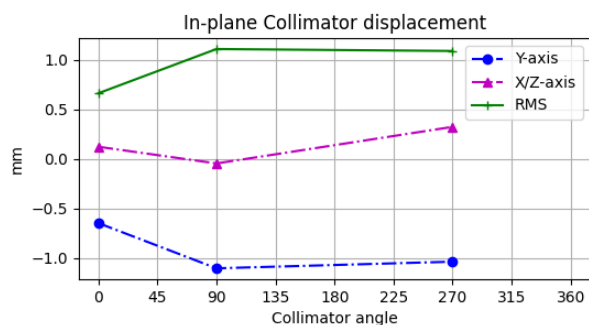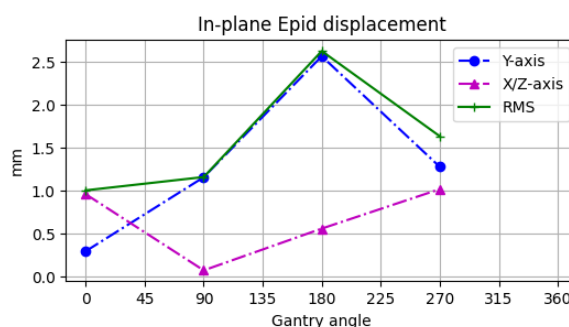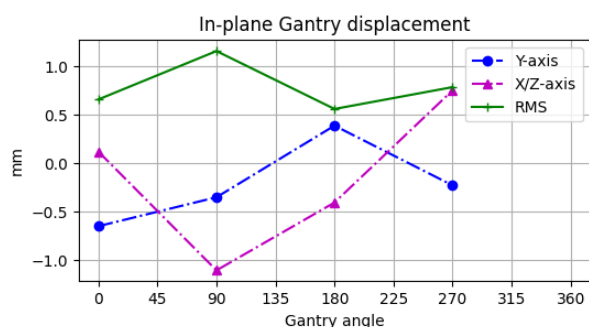

Notes:

Análise de Winston-Lutz com pylinac

Campo 2 cm X 2 cm

Tolerância: 1 mm

Limite de ação: 0.8 mm

Local das imagens: ./WL\WL_2019-09-27_analisado



Winston-Lutz Analysis

Metadata:


Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

Gantry images

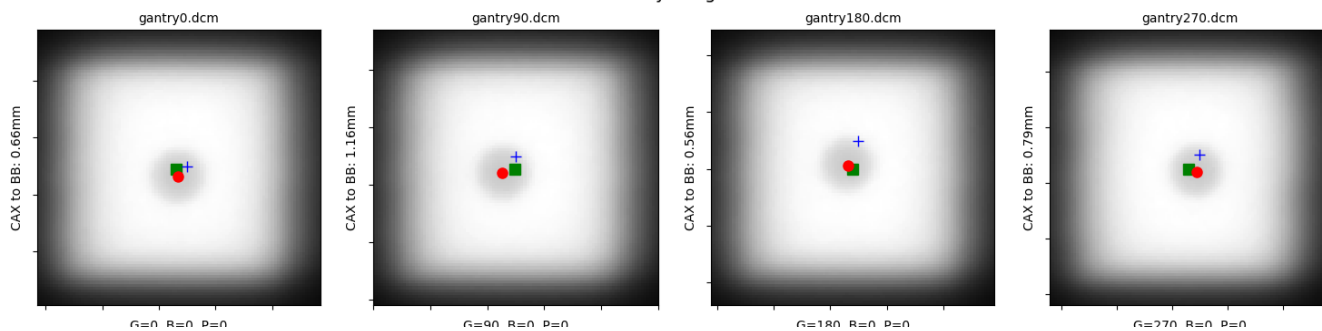



Winston-Lutz Analysis

Metadata:


Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

Collimator images

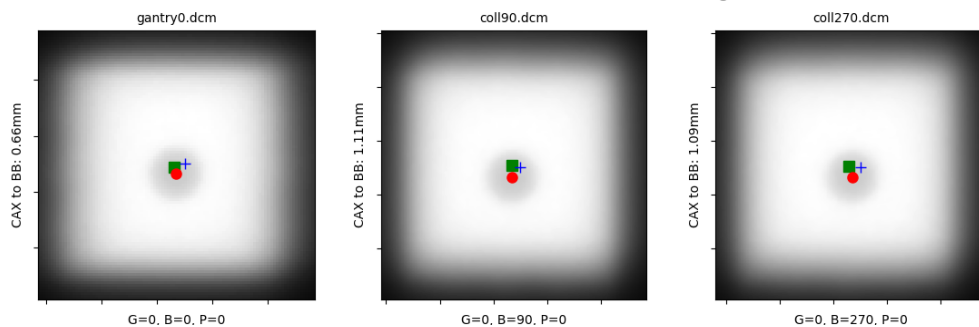



Winston-Lutz Analysis

Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

Couch images

