

Metadata: Winston-Lutz Analysis

Função: Físico Number of images: 8

Acelerador Linear: Synergy Full Maximum 2D CAX->BB distance: 1.76mm Modelo MLC: Agility Median 2D CAX->BB distance: 0.91mm

Shift to iso: facing gantry, move BB: LEFT 0.70mm; OUT 0.48mm; DOWN 0.67mm

Gantry 3D isocenter diameter: 1.38mm (4/8 images considered)

Maximum Gantry RMS deviation (mm): 1.57mm Maximum EPID RMS deviation (mm): 3.05mm


Gantry+Collimator 3D isocenter diameter: 1.40mm (6/8 images considered)

Collimator 2D isocenter diameter: 0.22mm (3/8 images considered)

Maximum Collimator RMS deviation (mm): 0.58

Couch 2D isocenter diameter: 2.50mm (3/8 images considered)

Maximum Couch RMS deviation (mm): 1.76

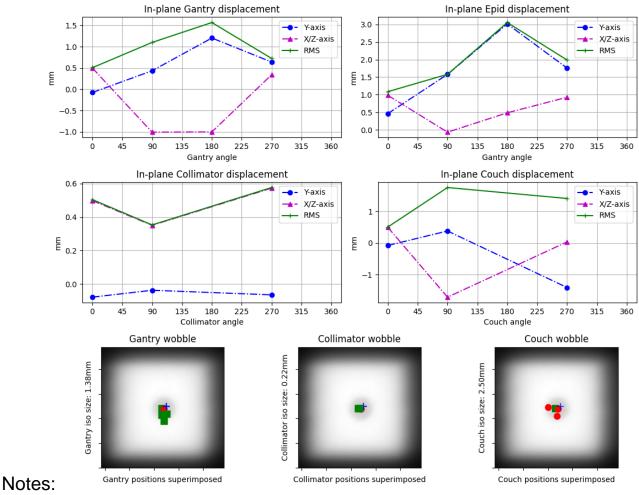

Análise de Winston-Lutz com pylinac

Campo 2 cm X 2 cm

Tolerância: 1 mm Limite de ação: 0.8 mm

Local das imagens: ./WL\WL\_2020-05-25\_analisado




Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

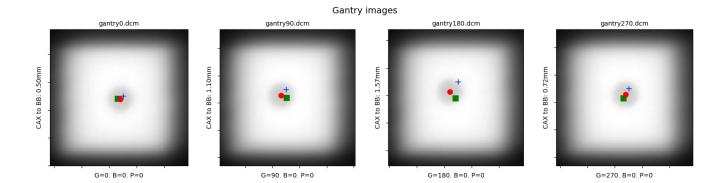




Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

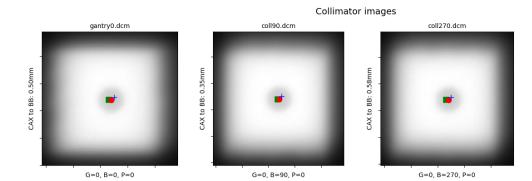



Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

