

Metadata: Winston-Lutz Analysis

Função: Físico Number of images: 8

Acelerador Linear: Synergy Full Maximum 2D CAX->BB distance: 0.91mm Modelo MLC: Agility Median 2D CAX->BB distance: 0.67mm

Shift to iso: facing gantry, move BB: RIGHT 0.36mm; IN 0.03mm; DOWN 0.24mm

Gantry 3D isocenter diameter: 1.14mm (4/8 images considered)

Maximum Gantry RMS deviation (mm): 0.74mm Maximum EPID RMS deviation (mm): 3.13mm


Gantry+Collimator 3D isocenter diameter: 1.41mm (6/8 images considered)

Collimator 2D isocenter diameter: 0.47mm (3/8 images considered)

Maximum Collimator RMS deviation (mm): 0.91

Couch 2D isocenter diameter: 0.80mm (3/8 images considered)

Maximum Couch RMS deviation (mm): 0.89

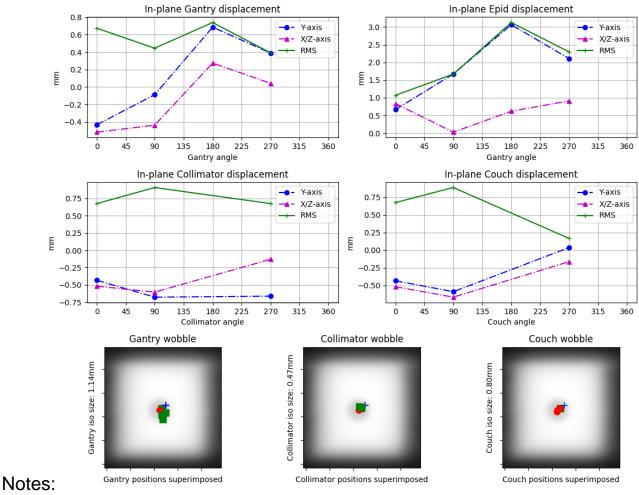

Análise de Winston-Lutz com pylinac

Campo 2 cm X 2 cm
Tolerância: 1 mm

Limite de ação: 0.8 mm

Local das imagens: ./WL\WL_2020-07-07_analisado




Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

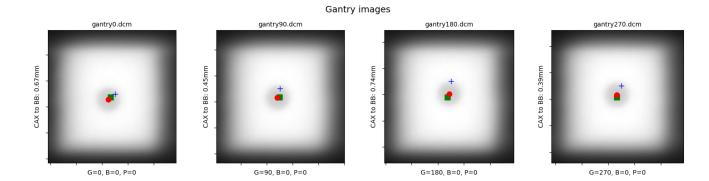




Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

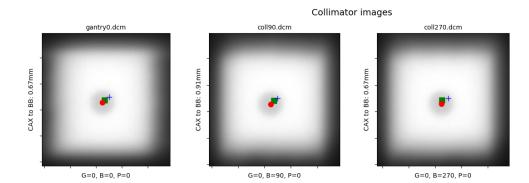



Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

