



# Winston-Lutz Analysis

## Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

## Winston-Lutz Analysis

Number of images: 8

Maximum 2D CAX->BB distance: 1.31mm

Median 2D CAX->BB distance: 0.92mm

Shift to iso: facing gantry, move BB: RIGHT 0.21mm; IN 0.19mm; DOWN 0.53mm

Gantry 3D isocenter diameter: 1.25mm (4/8 images considered)

Maximum Gantry RMS deviation (mm): 0.84mm

Maximum EPID RMS deviation (mm): 2.86mm








Gantry+Collimator 3D isocenter diameter: 1.33mm (6/8 images considered)

Collimator 2D isocenter diameter: 0.12mm (3/8 images considered)

Maximum Collimator RMS deviation (mm): 1.03

Couch 2D isocenter diameter: 1.28mm (3/8 images considered)

Maximum Couch RMS deviation (mm): 1.31

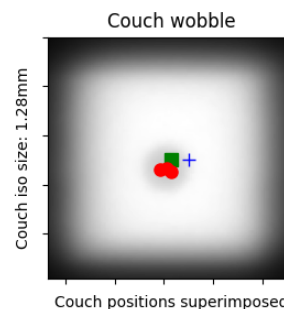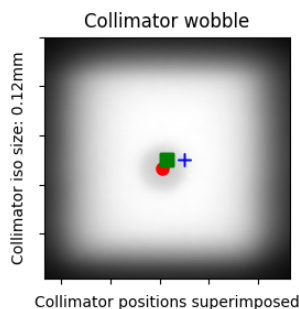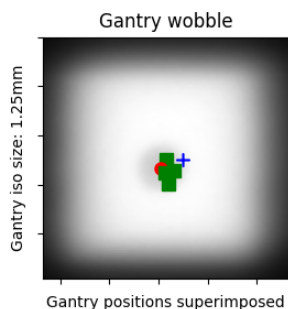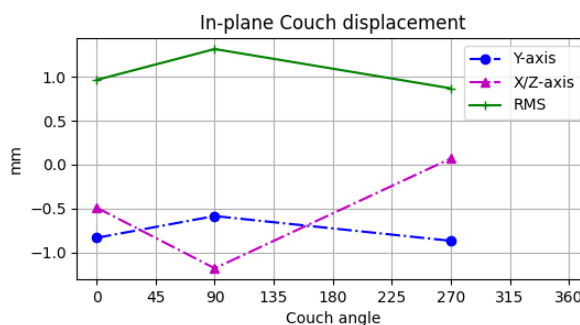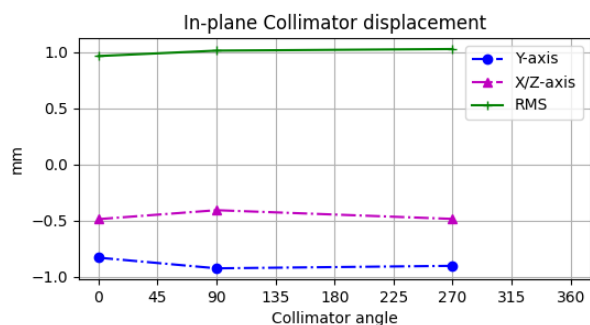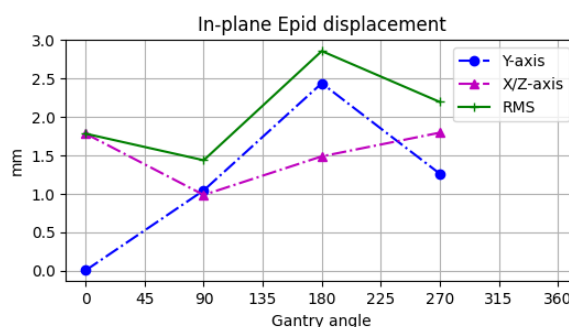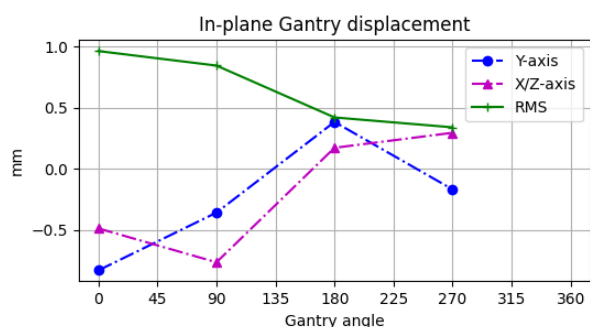

## Notes:

Análise de Winston-Lutz com pylinac

Campo 2 cm X 2 cm

Tolerância: 1 mm

Limite de ação: 0.8 mm

Local das imagens: ./WL\WL\_2019-11-06\_analisado



# Winston-Lutz Analysis

## Metadata:


Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

Gantry images

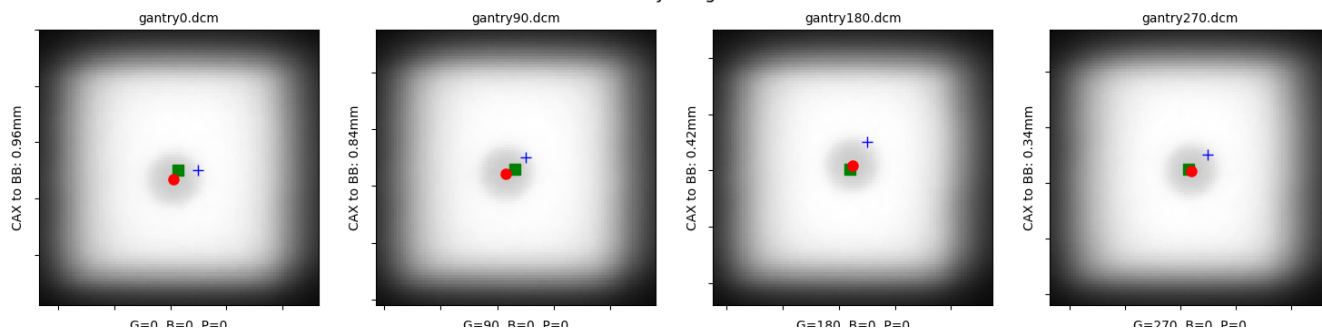



# Winston-Lutz Analysis

## Metadata:


Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

Collimator images

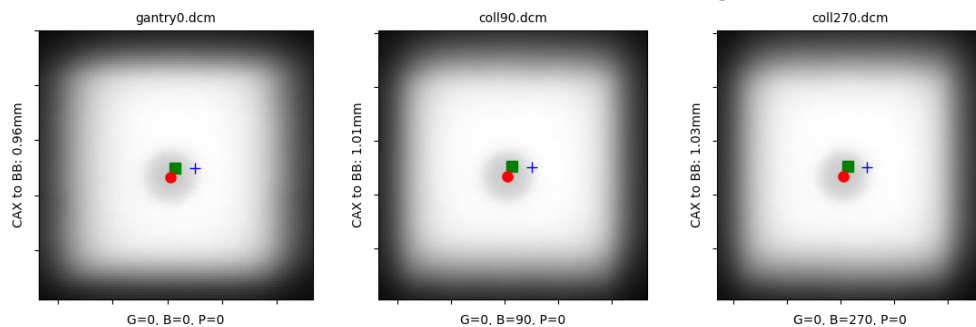



# Winston-Lutz Analysis

## Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

Couch images

