

Metadata: Winston-Lutz Analysis

Função: Físico Number of images: 8

Acelerador Linear: Synergy Full Maximum 2D CAX->BB distance: 1.45mm Modelo MLC: Agility Median 2D CAX->BB distance: 0.87mm

Shift to iso: facing gantry, move BB: RIGHT 0.56mm; IN 0.21mm; DOWN 0.53mm

Gantry 3D isocenter diameter: 1.07mm (4/8 images considered)

Maximum Gantry RMS deviation (mm): 0.75mm Maximum EPID RMS deviation (mm): 2.72mm


Gantry+Collimator 3D isocenter diameter: 1.62mm (6/8 images considered)

Collimator 2D isocenter diameter: 0.61mm (3/8 images considered)

Maximum Collimator RMS deviation (mm): 1.32

Couch 2D isocenter diameter: 1.03mm (3/8 images considered)

Maximum Couch RMS deviation (mm): 1.45

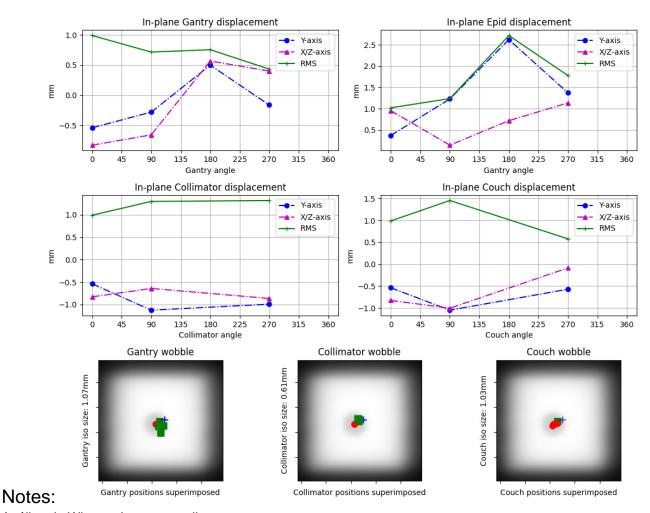

Análise de Winston-Lutz com pylinac

Campo 2 cm X 2 cm Tolerância: 1 mm Limite de ação: 0.8 mm

Local das imagens: ./WL\WL\_2019-09-04\_analisado




Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

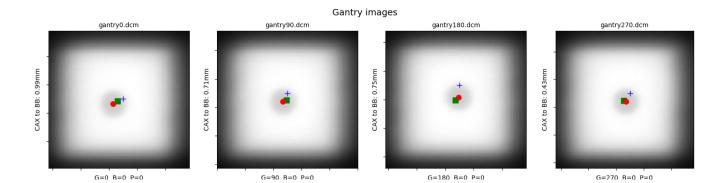




Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

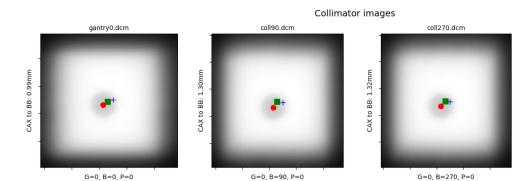



Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

