

Metadata: Winston-Lutz Analysis

Função: Físico Number of images: 9

Acelerador Linear: Synergy Full Maximum 2D CAX->BB distance: 1.15mm Modelo MLC: Agility Median 2D CAX->BB distance: 0.42mm

Shift to iso: facing gantry, move BB: RIGHT 0.04mm; OUT 0.30mm; DOWN 0.47mm

Gantry 3D isocenter diameter: 1.11mm (4/9 images considered)

Maximum Gantry RMS deviation (mm): 1.15mm Maximum EPID RMS deviation (mm): 34.12mm


Gantry+Collimator 3D isocenter diameter: 1.55mm (7/9 images considered)

Collimator 2D isocenter diameter: 0.44mm (4/9 images considered)

Maximum Collimator RMS deviation (mm): 0.42

Couch 2D isocenter diameter: 1.23mm (3/9 images considered)

Maximum Couch RMS deviation (mm): 0.92

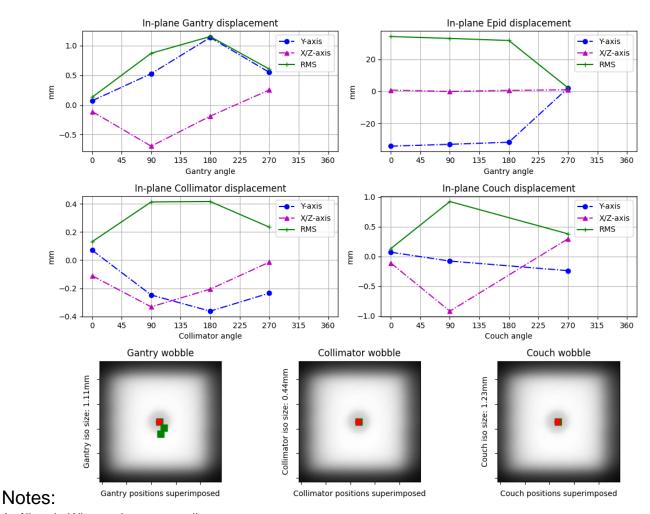

Análise de Winston-Lutz com pylinac

Campo 2 cm X 2 cm Tolerância: 1 mm Limite de ação: 0.8 mm

Local das imagens: ./WL_IOSP_20201216\WL_2020-11-18




Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

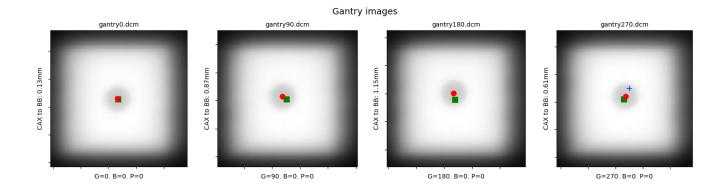




Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

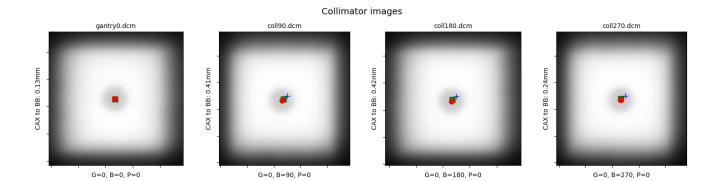



Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

