

Metadata: Winston-Lutz Analysis

Função: Físico Number of images: 8

Acelerador Linear: Synergy Full Maximum 2D CAX->BB distance: 1.67mm Modelo MLC: Agility Median 2D CAX->BB distance: 0.86mm

Shift to iso: facing gantry, move BB: LEFT 0.12mm; OUT 0.14mm; DOWN 1.24mm

Gantry 3D isocenter diameter: 1.06mm (4/8 images considered)

Maximum Gantry RMS deviation (mm): 1.54mm Maximum EPID RMS deviation (mm): 2.82mm


Gantry+Collimator 3D isocenter diameter: 1.40mm (6/8 images considered)

Collimator 2D isocenter diameter: 0.54mm (3/8 images considered)

Maximum Collimator RMS deviation (mm): 0.68

Couch 2D isocenter diameter: 1.81mm (3/8 images considered)

Maximum Couch RMS deviation (mm): 1.67

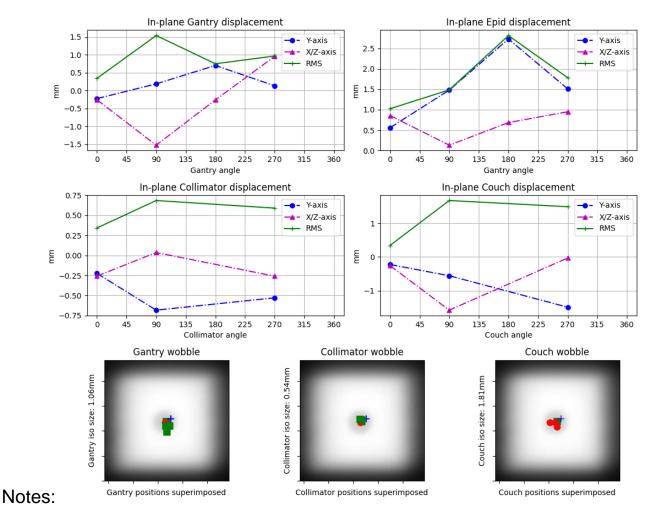

Análise de Winston-Lutz com pylinac

Campo 2 cm X 2 cm Tolerância: 1 mm Limite de ação: 0.8 mm

Local das imagens: ./WL\WL_2020-05-21_analisado




Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

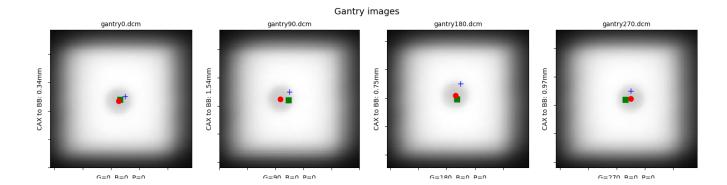




Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

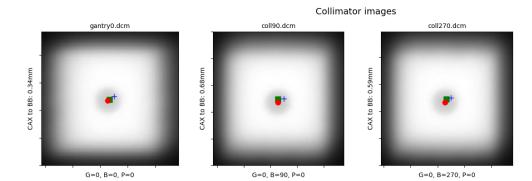



Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

