

Metadata: Winston-Lutz Analysis

Autor: William A. P. dos Santos

Função: Físico Number of images: 8

Acelerador Linear: Synergy Full Maximum 2D CAX->BB distance: 1.21mm Modelo MLC: Agility Median 2D CAX->BB distance: 0.84mm

Shift to iso: facing gantry, move BB: RIGHT 0.41mm; IN 0.28mm; DOWN 0.38mm

Gantry 3D isocenter diameter: 1.15mm (4/8 images considered)

Maximum Gantry RMS deviation (mm): 0.69mm Maximum EPID RMS deviation (mm): 2.61mm


Gantry+Collimator 3D isocenter diameter: 1.43mm (6/8 images considered)

Collimator 2D isocenter diameter: 0.34mm (3/8 images considered)

Maximum Collimator RMS deviation (mm): 1.21

Couch 2D isocenter diameter: 0.80mm (3/8 images considered)

Maximum Couch RMS deviation (mm): 0.99

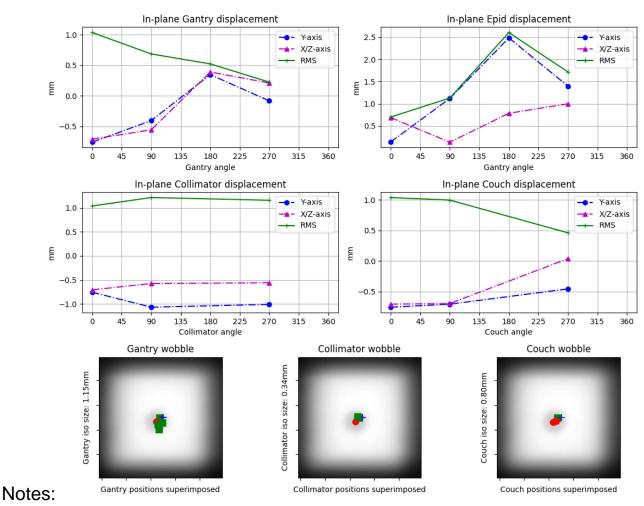

Análise de Winston-Lutz com pylinac

Campo 2 cm X 2 cm Tolerância: 1 mm Limite de ação: 0.8 mm

Local das imagens: ./WL\WL_2020-08-05_analisado




Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

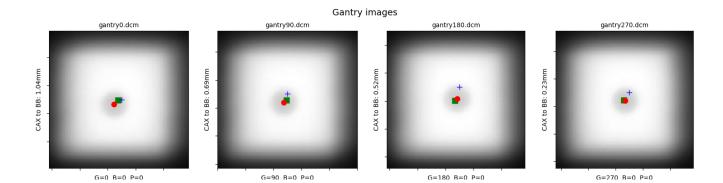




Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

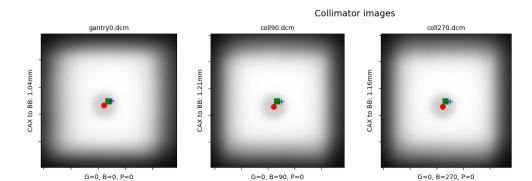



Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

