



# Winston-Lutz Analysis

## Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

## Winston-Lutz Analysis

Number of images: 9

Maximum 2D CAX->BB distance: 0.98mm

Median 2D CAX->BB distance: 0.74mm

Shift to iso: facing gantry, move BB: RIGHT 0.20mm; OUT 0.11mm; DOWN 0.52mm

Gantry 3D isocenter diameter: 1.47mm (4/9 images considered)

Maximum Gantry RMS deviation (mm): 0.97mm

Maximum EPID RMS deviation (mm): 2.97mm








Gantry+Collimator 3D isocenter diameter: 1.52mm (7/9 images considered)

Collimator 2D isocenter diameter: 0.49mm (4/9 images considered)

Maximum Collimator RMS deviation (mm): 0.74

Couch 2D isocenter diameter: 0.94mm (3/9 images considered)

Maximum Couch RMS deviation (mm): 0.98

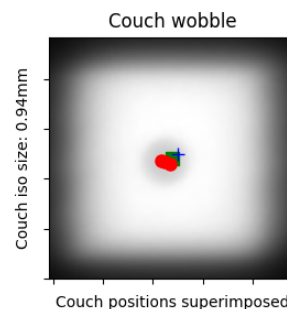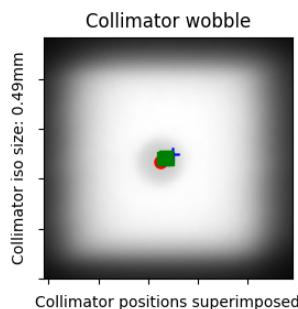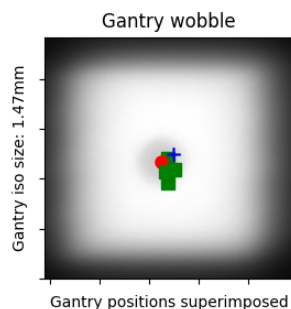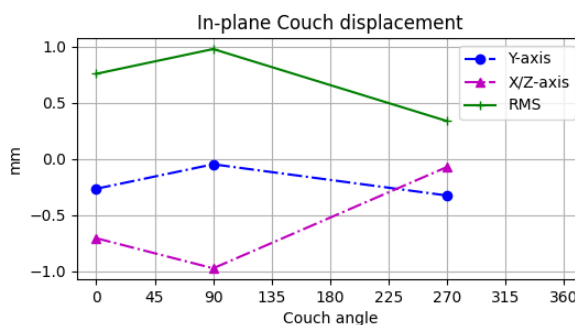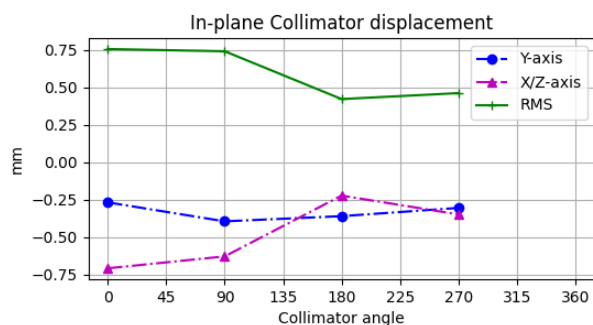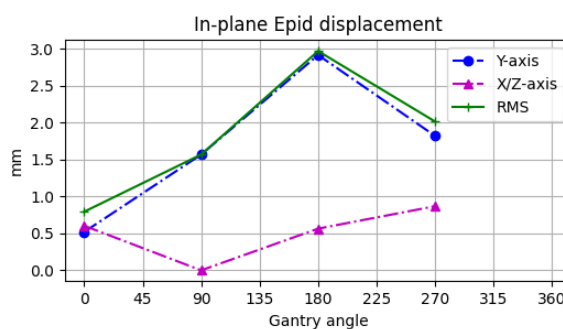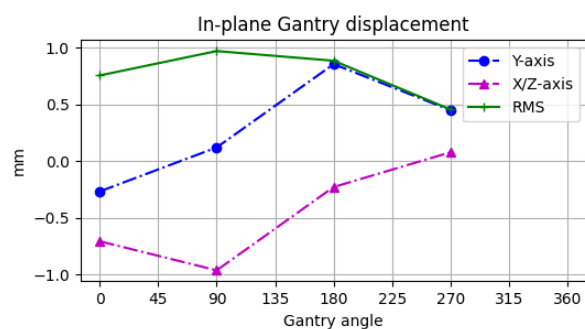

## Notes:

Análise de Winston-Lutz com pylinac

Campo 2 cm X 2 cm

Tolerância: 1 mm

Limite de ação: 0.8 mm

Local das imagens: ./WL\_IOSP\_20201216\WL\_2020-10-26



# Winston-Lutz Analysis

## Metadata:


Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

Gantry images

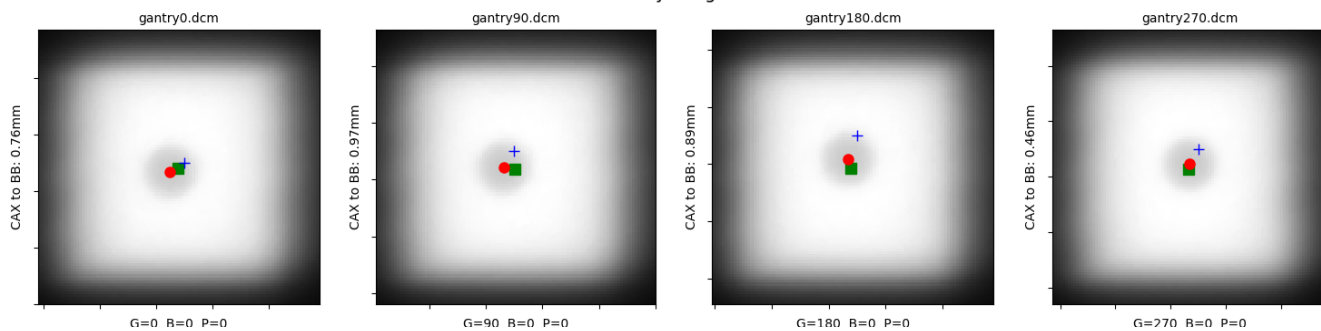



# Winston-Lutz Analysis

## Metadata:


Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

Collimator images

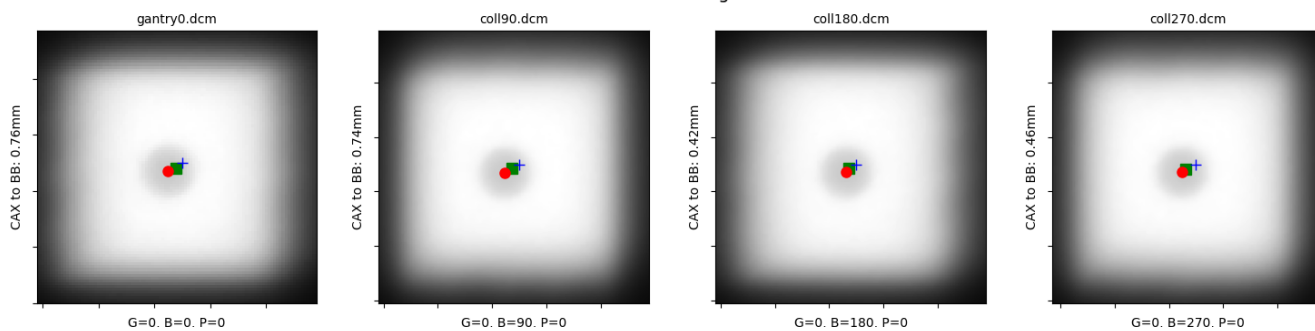



# Winston-Lutz Analysis

## Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

Couch images

