



Winston-Lutz Analysis

Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

Winston-Lutz Analysis

Number of images: 9

Maximum 2D CAX->BB distance: 1.13mm

Median 2D CAX->BB distance: 0.77mm

Shift to iso: facing gantry, move BB: LEFT 0.20mm; OUT 0.03mm; UP 0.37mm

Gantry 3D isocenter diameter: 1.22mm (4/9 images considered)

Maximum Gantry RMS deviation (mm): 0.99mm

Maximum EPID RMS deviation (mm): 2.75mm








Gantry+Collimator 3D isocenter diameter: 1.75mm (7/9 images considered)

Collimator 2D isocenter diameter: 0.54mm (4/9 images considered)

Maximum Collimator RMS deviation (mm): 0.82

Couch 2D isocenter diameter: 1.55mm (3/9 images considered)

Maximum Couch RMS deviation (mm): 1.13

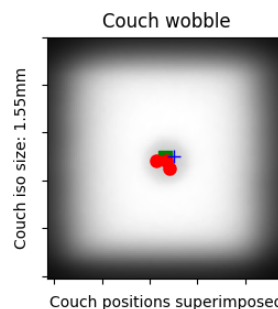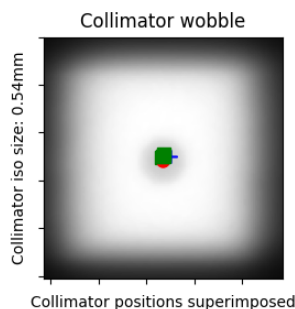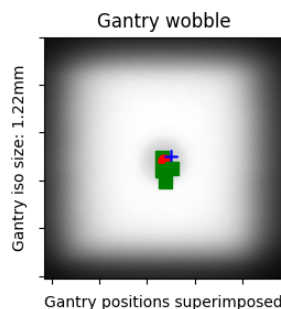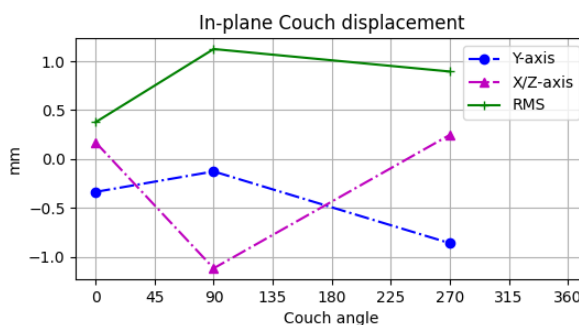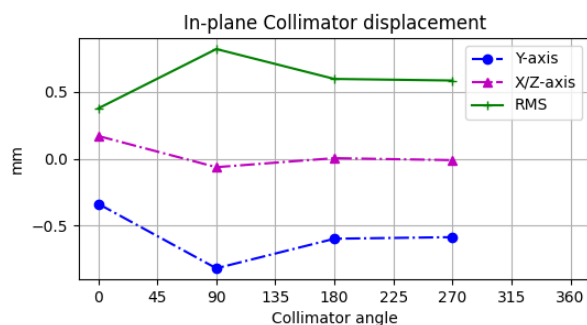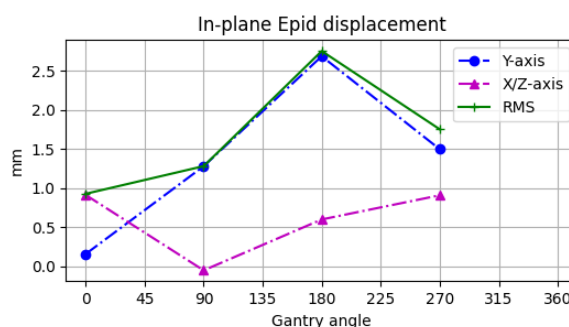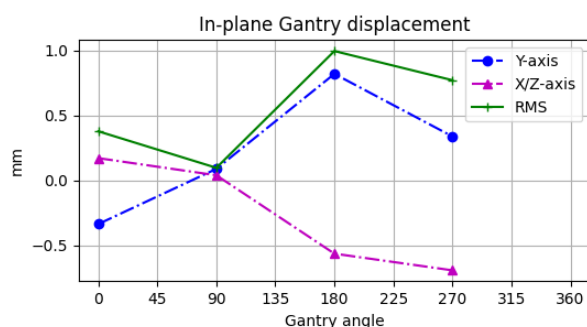

Notes:

Análise de Winston-Lutz com pylinac

Campo 2 cm X 2 cm

Tolerância: 1 mm

Limite de ação: 0.8 mm

Local das imagens: ./WL_IOSP_20201216\WL_2020-12-11



Winston-Lutz Analysis

Metadata:


Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

Gantry images

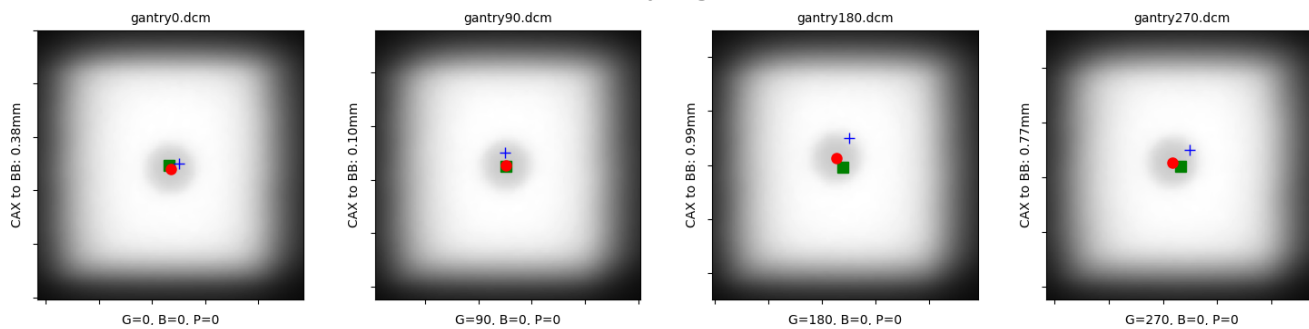



Winston-Lutz Analysis

Metadata:


Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

Collimator images

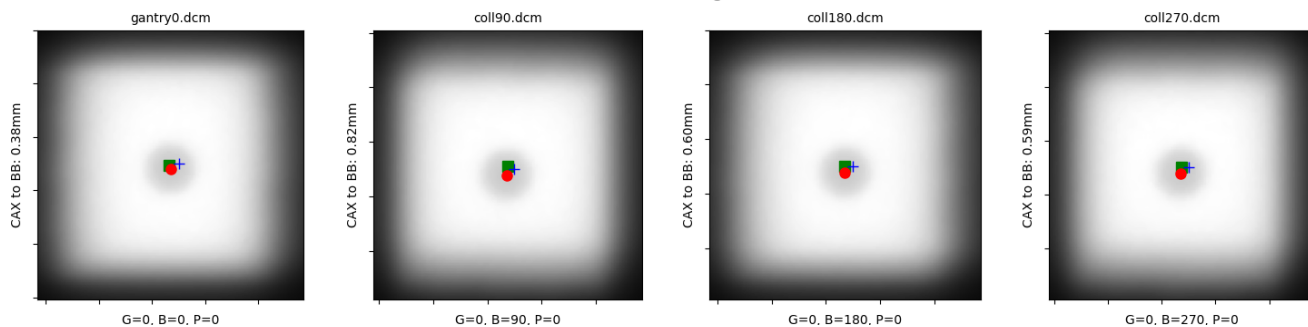



Winston-Lutz Analysis

Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

Couch images

