

Metadata: Winston-Lutz Analysis

Função: Físico Number of images: 8

Acelerador Linear: Synergy Full Maximum 2D CAX->BB distance: 1.36mm Modelo MLC: Agility Median 2D CAX->BB distance: 0.65mm

Shift to iso: facing gantry, move BB: RIGHT 0.03mm; IN 0.02mm; DOWN 0.27mm

Gantry 3D isocenter diameter: 1.31mm (4/8 images considered)

Maximum Gantry RMS deviation (mm): 0.67mm Maximum EPID RMS deviation (mm): 1.70mm


Gantry+Collimator 3D isocenter diameter: 1.32mm (6/8 images considered)

Collimator 2D isocenter diameter: 0.12mm (3/8 images considered)

Maximum Collimator RMS deviation (mm): 0.65

Couch 2D isocenter diameter: 1.11mm (3/8 images considered)

Maximum Couch RMS deviation (mm): 1.36

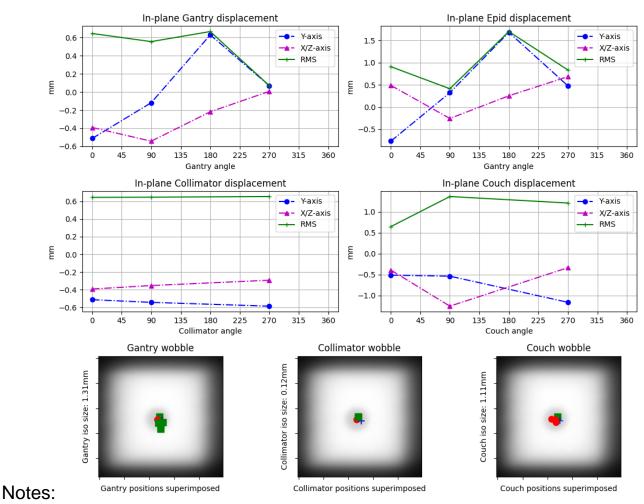

Análise de Winston-Lutz com pylinac

Campo 2 cm X 2 cm Tolerância: 1 mm Limite de ação: 0.8 mm

Local das imagens: ./WL\WL\_2020-06-15\_analisado




Metadata:


Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

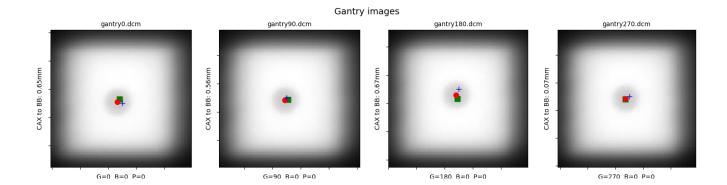

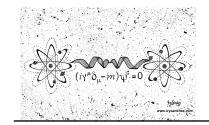


Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

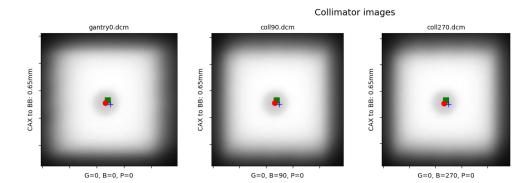



Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

