

Metadata: Winston-Lutz Analysis

Função: Físico Number of images: 9

Acelerador Linear: Synergy Full Maximum 2D CAX->BB distance: 1.22mm Modelo MLC: Agility Median 2D CAX->BB distance: 0.53mm

Shift to iso: facing gantry, move BB: RIGHT 0.28mm; OUT 0.34mm; DOWN 0.16mm

Gantry 3D isocenter diameter: 1.17mm (4/9 images considered)

Maximum Gantry RMS deviation (mm): 1.22mm Maximum EPID RMS deviation (mm): 3.08mm


Gantry+Collimator 3D isocenter diameter: 1.51mm (7/9 images considered)

Collimator 2D isocenter diameter: 0.37mm (4/9 images considered)

Maximum Collimator RMS deviation (mm): 0.53

Couch 2D isocenter diameter: 1.38mm (3/9 images considered)

Maximum Couch RMS deviation (mm): 1.05

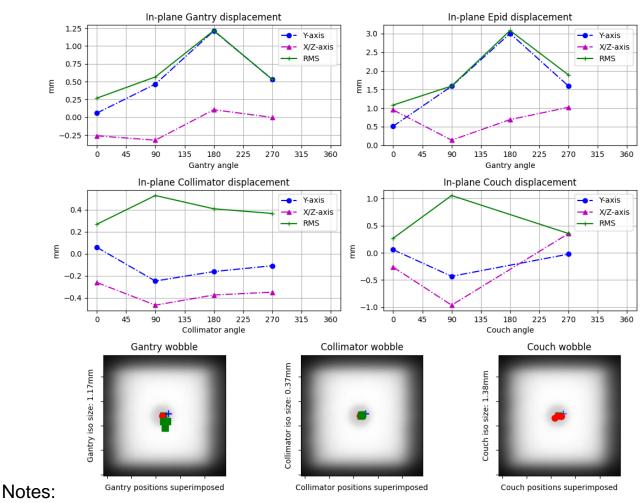

Análise de Winston-Lutz com pylinac

Campo 2 cm X 2 cm Tolerância: 1 mm Limite de ação: 0.8 mm

Local das imagens: ./WL_IOSP_20201216\WL_2020-11-13




Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

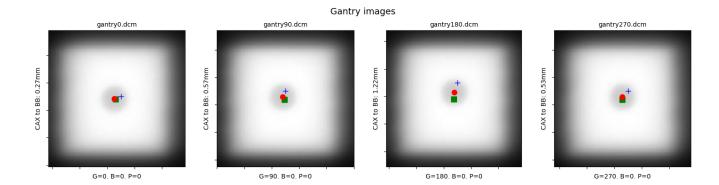




Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

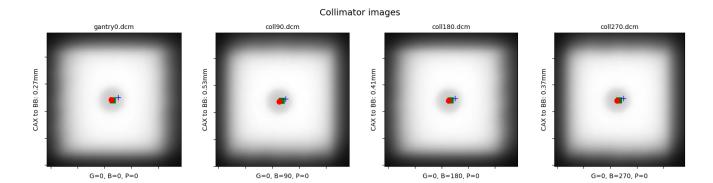



Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

