

Winston-Lutz Analysis

Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

Winston-Lutz Analysis

Number of images: 8

Maximum 2D CAX->BB distance: 1.26mm

Median 2D CAX->BB distance: 0.87mm

Shift to iso: facing gantry, move BB: RIGHT 0.18mm; IN 0.24mm; DOWN 0.47mm

Gantry 3D isocenter diameter: 1.20mm (4/8 images considered)

Maximum Gantry RMS deviation (mm): 0.82mm

Maximum EPID RMS deviation (mm): 3.91mm






Gantry+Collimator 3D isocenter diameter: 1.62mm (6/8 images considered)

Collimator 2D isocenter diameter: 0.39mm (3/8 images considered)


Maximum Collimator RMS deviation (mm): 1.23

Couch 2D isocenter diameter: 1.18mm (3/8 images considered)


Maximum Couch RMS deviation (mm): 1.26

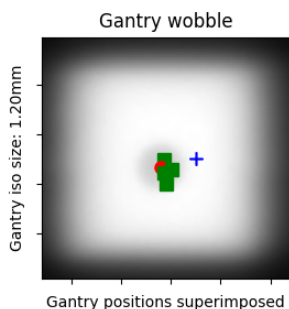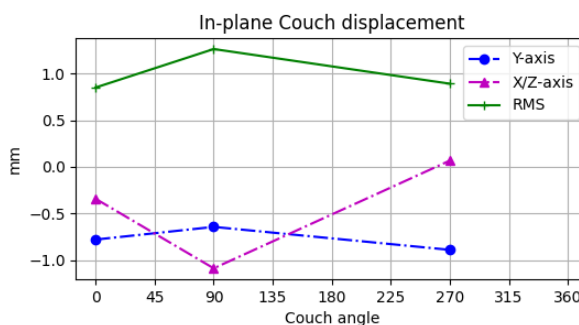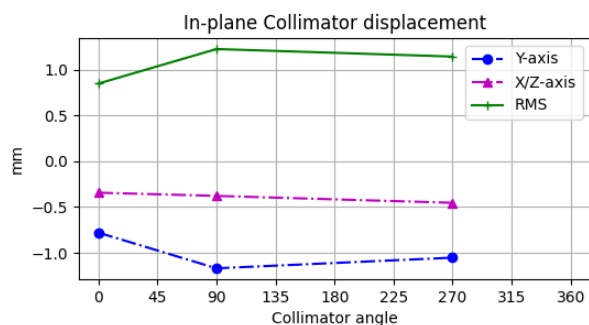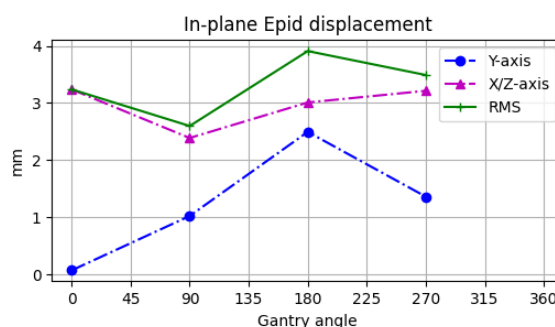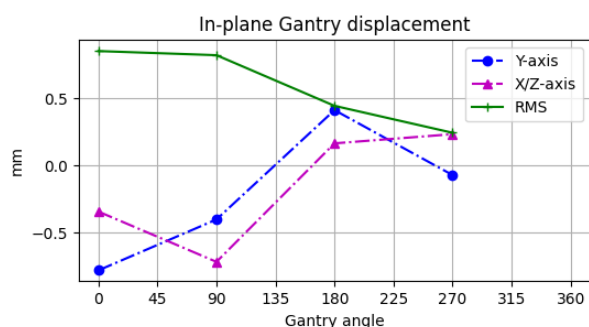

Gantry positions superimposed

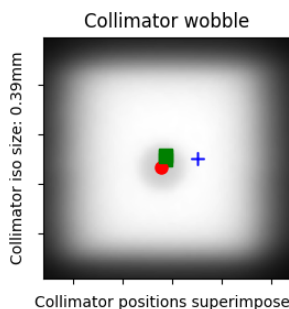

Collimator positions superimposed

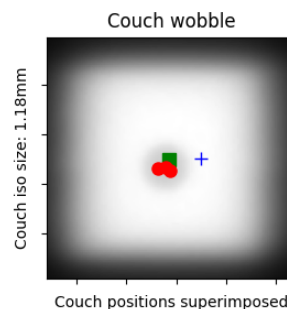

Couch positions superimposed

Notes:

Análise de Winston-Lutz com pylinac

Campo 2 cm X 2 cm

Tolerância: 1 mm

Limite de ação: 0.8 mm

Local das imagens: ./WL\WL_2019-11-12_analisado



Winston-Lutz Analysis

Metadata:


Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

Gantry images

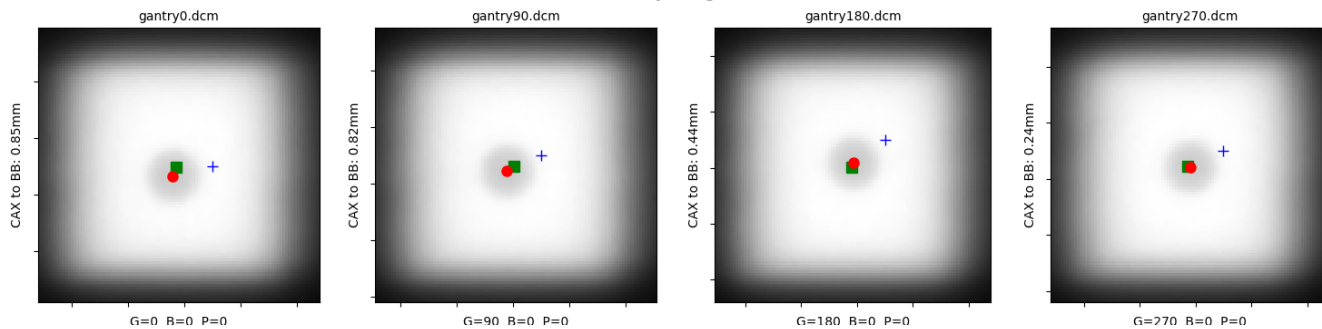



Winston-Lutz Analysis

Metadata:


Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

Collimator images

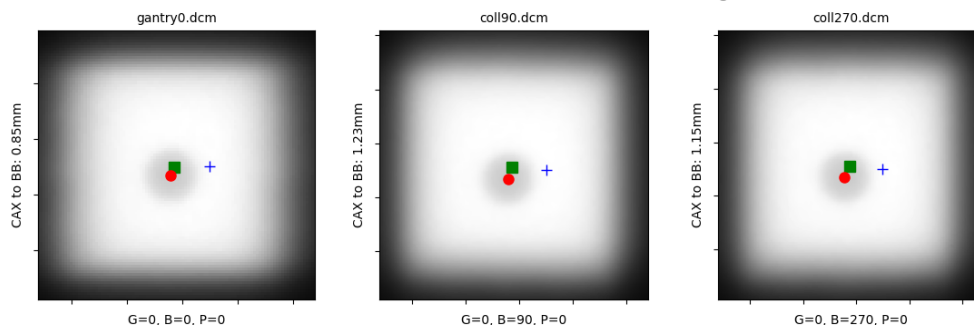



Winston-Lutz Analysis

Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

Couch images

