

Metadata: Winston-Lutz Analysis

Autor: William A. P. dos Santos

Função: Físico Number of images: 8

Acelerador Linear: Synergy Full Maximum 2D CAX->BB distance: 1.13mm Modelo MLC: Agility Median 2D CAX->BB distance: 0.78mm

Shift to iso: facing gantry, move BB: RIGHT 0.30mm; OUT 0.03mm; DOWN 0.94mm

Gantry 3D isocenter diameter: 1.29mm (4/8 images considered)

Maximum Gantry RMS deviation (mm): 1.13mm Maximum EPID RMS deviation (mm): 2.90mm


Gantry+Collimator 3D isocenter diameter: 1.52mm (6/8 images considered)

Collimator 2D isocenter diameter: 0.28mm (3/8 images considered)

Maximum Collimator RMS deviation (mm): 0.86

Couch 2D isocenter diameter: 1.00mm (3/8 images considered)

Maximum Couch RMS deviation (mm): 1.10

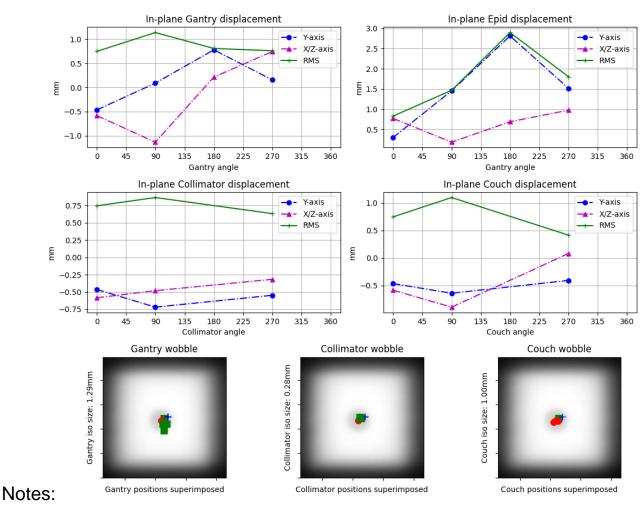

Análise de Winston-Lutz com pylinac

Campo 2 cm X 2 cm Tolerância: 1 mm Limite de ação: 0.8 mm

Local das imagens: ./WL\WL_2019-07-29_analisado




Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

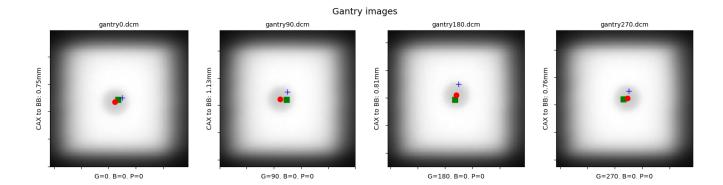




Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

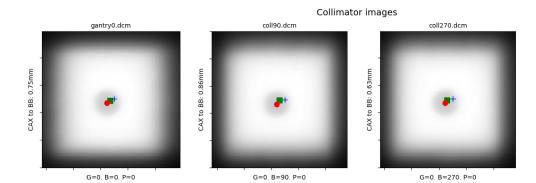



Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

