



# Winston-Lutz Analysis

## Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

## Winston-Lutz Analysis

Number of images: 9

Maximum 2D CAX->BB distance: 0.93mm

Median 2D CAX->BB distance: 0.47mm

Shift to iso: facing gantry, move BB: LEFT 0.12mm; OUT 0.17mm; DOWN 0.16mm

Gantry 3D isocenter diameter: 1.20mm (4/9 images considered)

Maximum Gantry RMS deviation (mm): 0.87mm

Maximum EPID RMS deviation (mm): 2.58mm






Gantry+Collimator 3D isocenter diameter: 1.20mm (7/9 images considered)

Collimator 2D isocenter diameter: 0.35mm (4/9 images considered)


Maximum Collimator RMS deviation (mm): 0.37

Couch 2D isocenter diameter: 1.29mm (3/9 images considered)


Maximum Couch RMS deviation (mm): 0.93

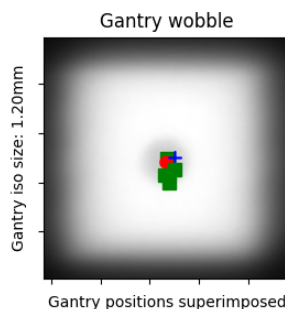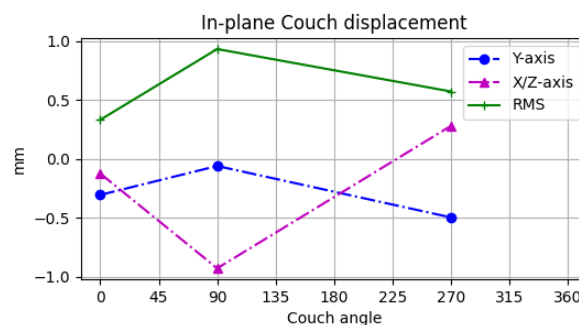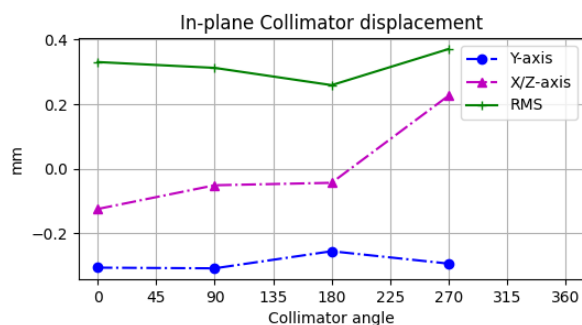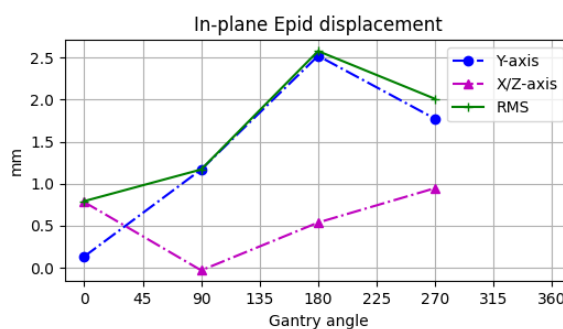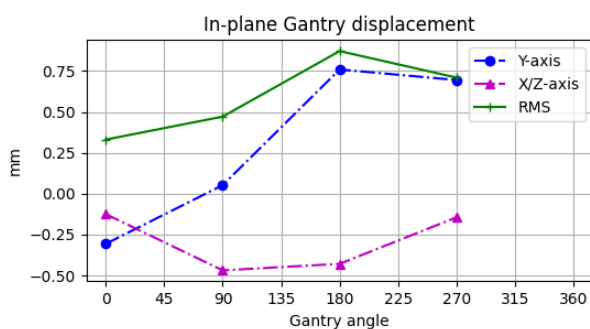

Gantry positions superimposed

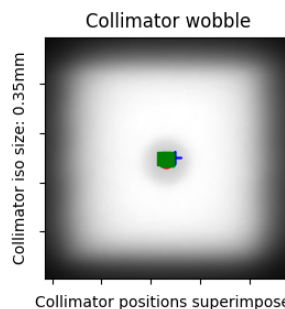

Collimator positions superimposed

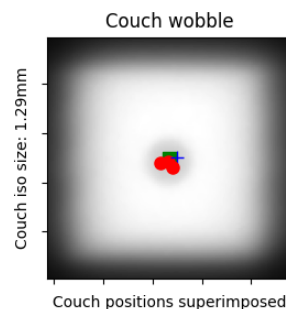

Couch positions superimposed

## Notes:

Análise de Winston-Lutz com pylinac

Campo 2 cm X 2 cm

Tolerância: 1 mm

Limite de ação: 0.8 mm

Local das imagens: ./WL\_IOSP\_20201216\WL\_2020-11-05



# Winston-Lutz Analysis

## Metadata:


Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

Gantry images

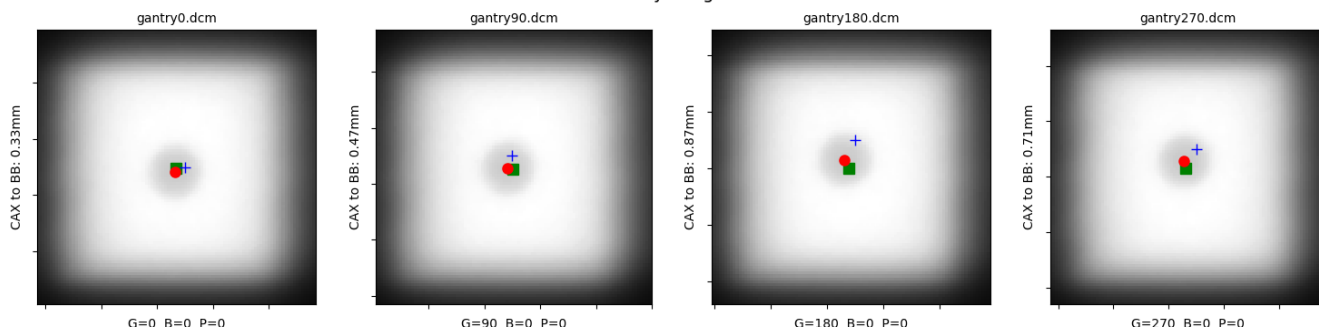



# Winston-Lutz Analysis

## Metadata:


Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

Collimator images

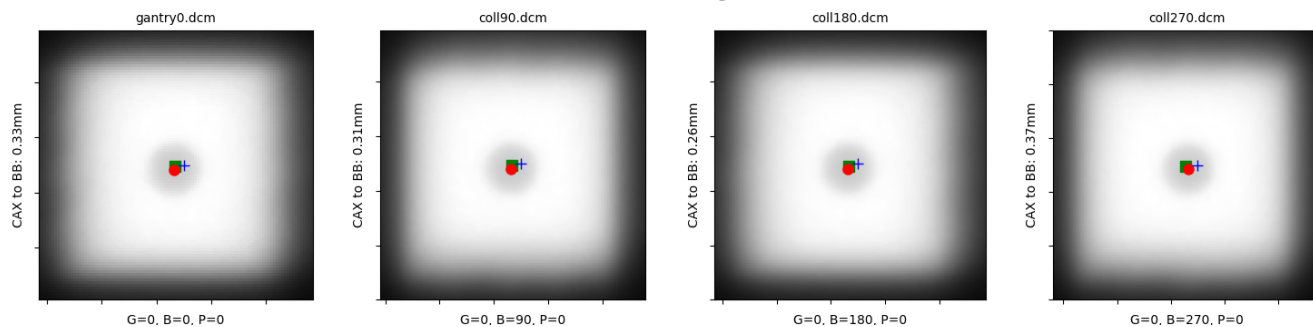



# Winston-Lutz Analysis

## Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

Couch images

