

Winston-Lutz Analysis

Metadata:

Winston-Lutz Analysis

Autor: William A. P. dos Santos =====

Função: Físico

Number of images: 63

Acelerador Linear: Synergy Full Maximum 2D CAX->BB distance: 1.28mm

Modelo MLC: Agility Median 2D CAX->BB distance: 0.57mm

Shift to iso: facing gantry, move BB: RIGHT 0.14mm; IN 0.24mm; DOWN 0.06mm

Gantry 3D isocenter diameter: 1.20mm (27/63 images considered)

Maximum Gantry RMS deviation (mm): 0.70mm

Maximum EPID RMS deviation (mm): 4.80mm


Gantry+Collimator 3D isocenter diameter: 1.34mm (51/63 images considered)

Collimator 2D isocenter diameter: 0.45mm (27/63 images considered)

Maximum Collimator RMS deviation (mm): 0.86

Couch 2D isocenter diameter: 1.58mm (15/63 images considered)

Maximum Couch RMS deviation (mm): 1.28

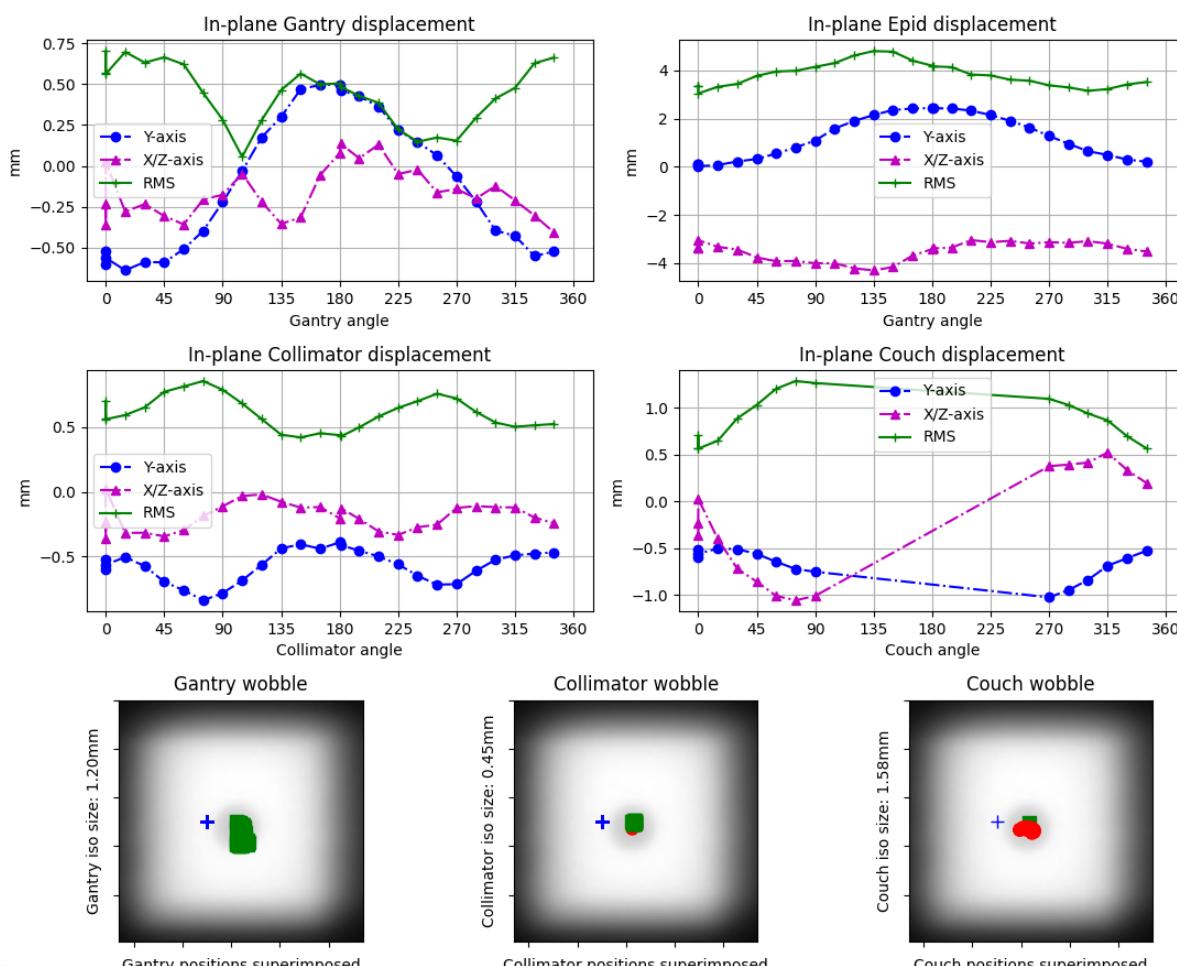

Notes:


Análise de Winston-Lutz com pylinac

Campo 2 cm X 2 cm

Tolerância: 1 mm

Límite de ação: 0.8 mm

Local das imagens: ./Teste_26AGO2020/imagem.jpg

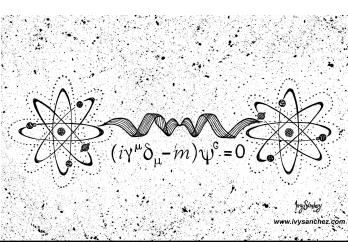

Winston-Lutz Analysis


Metadata:


Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synei

Modelo MLC: Agility

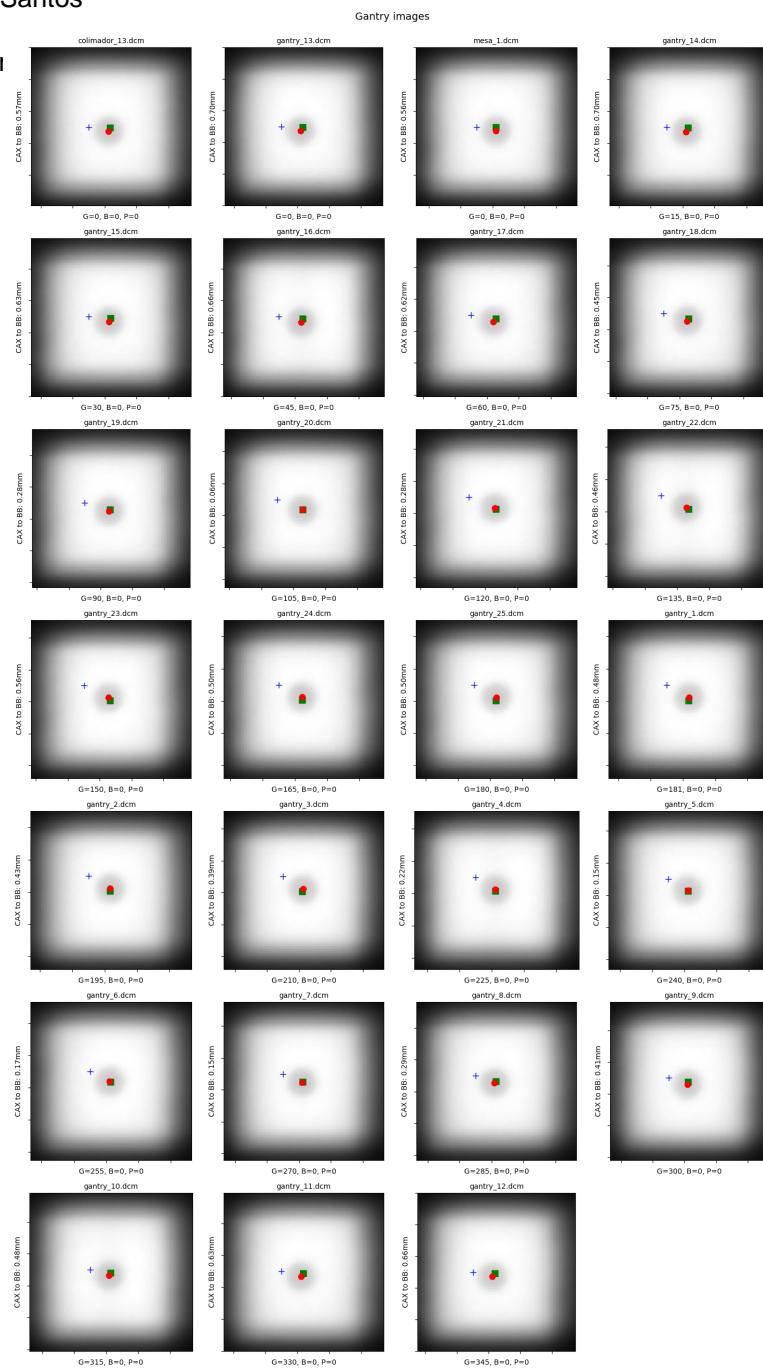

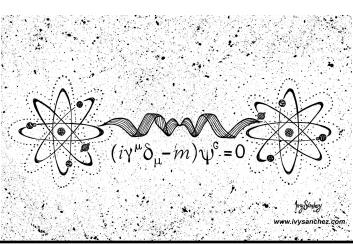

Winston-Lutz Analysis


Metadata:


Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synei

Modelo MLC: Agility

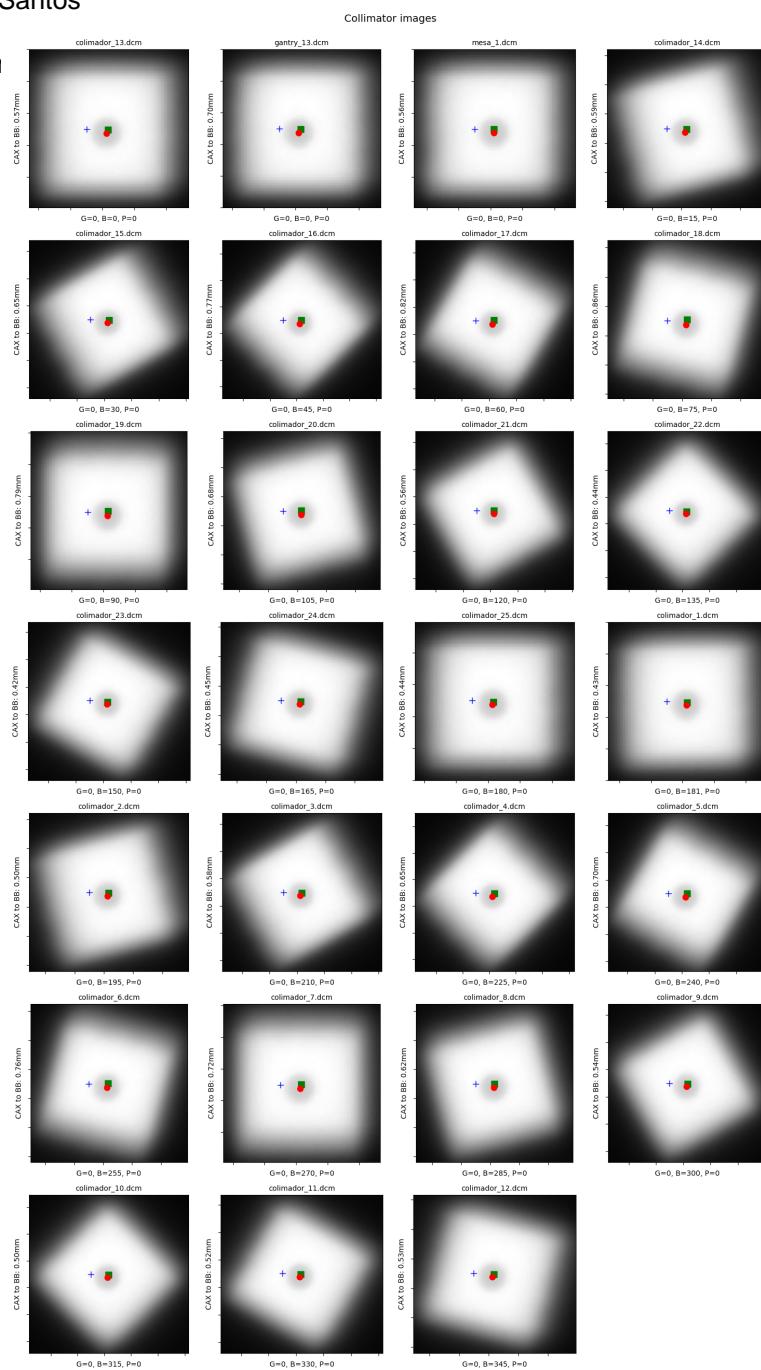

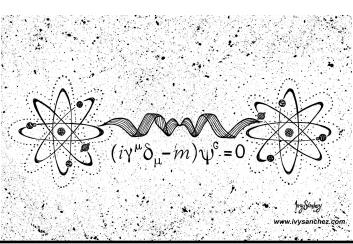

Winston-Lutz Analysis

Metadata:

Autor: William A. P. dos Santos

Couch images

