

Metadata: Winston-Lutz Analysis

Função: Físico Number of images: 7

Acelerador Linear: Synergy Full Maximum 2D CAX->BB distance: 1.04mm Modelo MLC: Agility Median 2D CAX->BB distance: 0.62mm

Shift to iso: facing gantry, move BB: RIGHT 0.41mm; IN 0.48mm; DOWN 0.09mm

Gantry 3D isocenter diameter: 1.05mm (4/7 images considered)

Maximum Gantry RMS deviation (mm): 0.44mm Maximum EPID RMS deviation (mm): 2.30mm


Gantry+Collimator 3D isocenter diameter: 1.20mm (6/7 images considered)

Collimator 2D isocenter diameter: 0.30mm (3/7 images considered)

Maximum Collimator RMS deviation (mm): 1.04

Couch 2D isocenter diameter: 0.57mm (2/7 images considered)

Maximum Couch RMS deviation (mm): 0.00

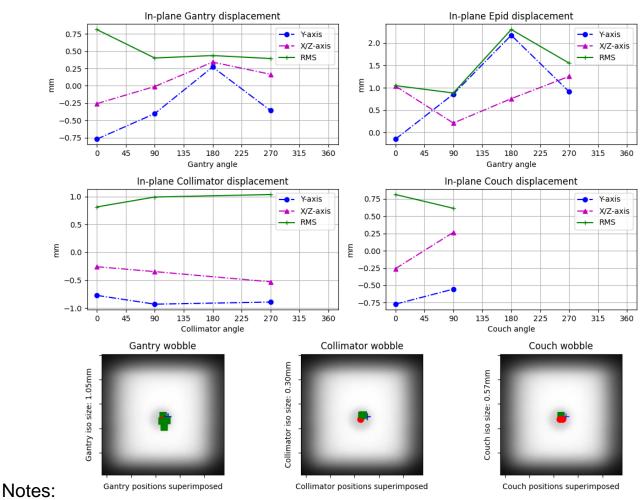

Análise de Winston-Lutz com pylinac

Campo 2 cm X 2 cm Tolerância: 1 mm Limite de ação: 0.8 mm

Local das imagens: ./WL\WL_2020-07-08_r




Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

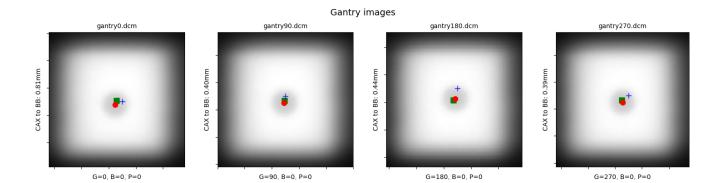




Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

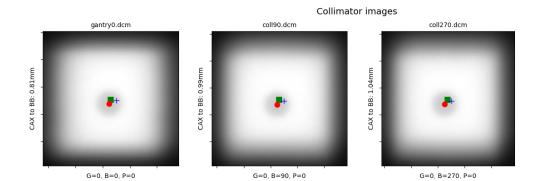



Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

