

Metadata: Winston-Lutz Analysis

Autor: William A. P. dos Santos

Função: Físico Number of images: 8

Acelerador Linear: Synergy Full Maximum 2D CAX->BB distance: 1.31mm Modelo MLC: Agility Median 2D CAX->BB distance: 0.61mm

Shift to iso: facing gantry, move BB: LEFT 0.24mm; OUT 0.01mm; DOWN 0.30mm

Gantry 3D isocenter diameter: 1.01mm (4/8 images considered)

Maximum Gantry RMS deviation (mm): 0.61mm Maximum EPID RMS deviation (mm): 2.62mm


Gantry+Collimator 3D isocenter diameter: 1.34mm (6/8 images considered)

Collimator 2D isocenter diameter: 0.40mm (3/8 images considered)

Maximum Collimator RMS deviation (mm): 0.81

Couch 2D isocenter diameter: 1.77mm (3/8 images considered)

Maximum Couch RMS deviation (mm): 1.31

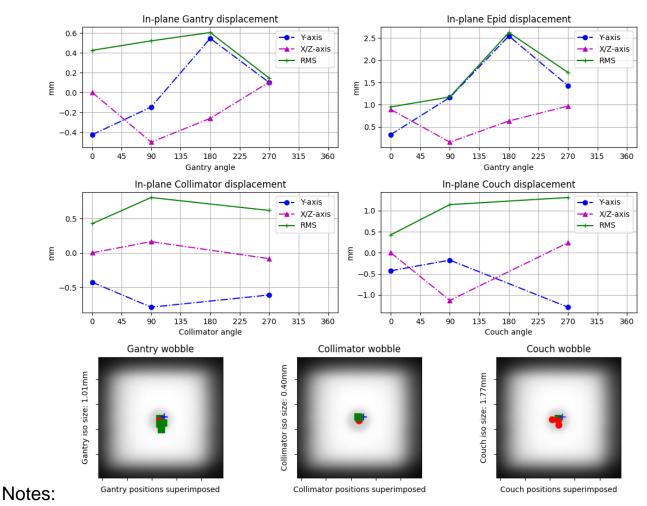

Análise de Winston-Lutz com pylinac

Campo 2 cm X 2 cm Tolerância: 1 mm Limite de ação: 0.8 mm

Local das imagens: ./WL\WL_2019-10-16_analisado




Metadata:


Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

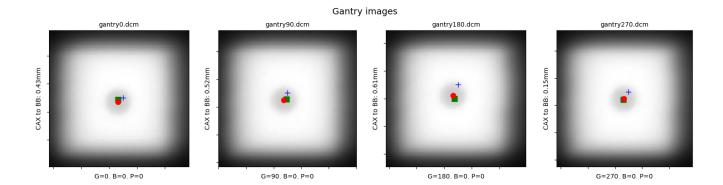

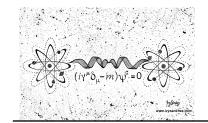


Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

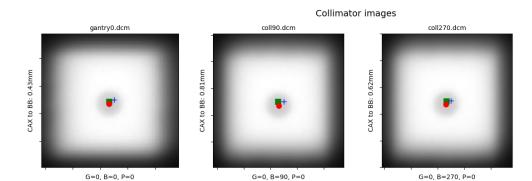



Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

