

Metadata: Winston-Lutz Analysis

Função: Físico Number of images: 8

Acelerador Linear: Synergy Full Maximum 2D CAX->BB distance: 1.40mm Modelo MLC: Agility Median 2D CAX->BB distance: 1.11mm

Shift to iso: facing gantry, move BB: RIGHT 0.11mm; IN 0.42mm; DOWN 0.27mm

Gantry 3D isocenter diameter: 1.06mm (4/8 images considered)

Maximum Gantry RMS deviation (mm): 0.79mm Maximum EPID RMS deviation (mm): 2.31mm


Gantry+Collimator 3D isocenter diameter: 1.51mm (6/8 images considered)

Collimator 2D isocenter diameter: 0.50mm (3/8 images considered)

Maximum Collimator RMS deviation (mm): 1.38

Couch 2D isocenter diameter: 1.07mm (3/8 images considered)

Maximum Couch RMS deviation (mm): 1.40

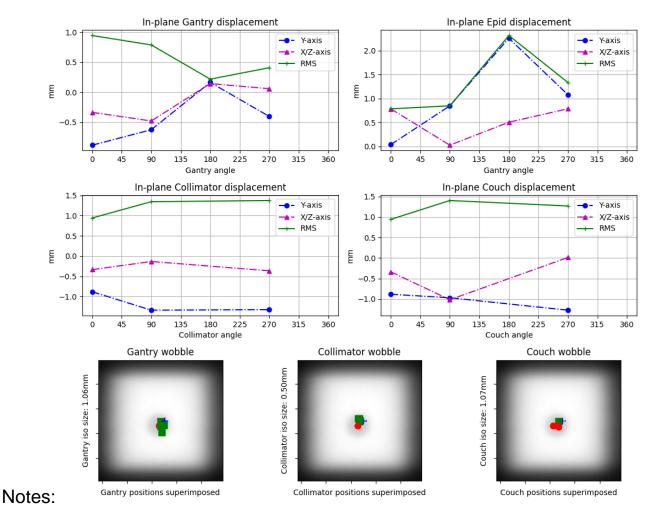

Análise de Winston-Lutz com pylinac

Campo 2 cm X 2 cm Tolerância: 1 mm Limite de ação: 0.8 mm

Local das imagens: ./WL\WL_2019-11-25_analisado




Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

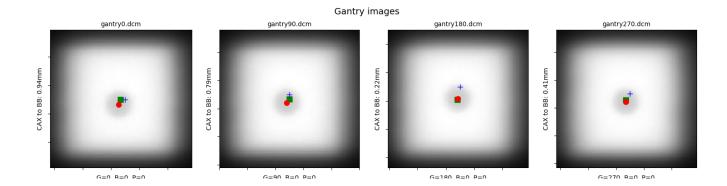




Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

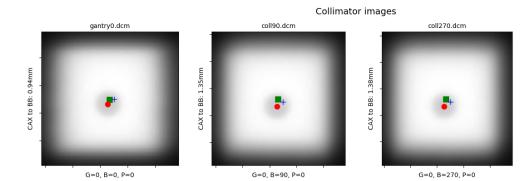



Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

