



Winston-Lutz Analysis

Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

Winston-Lutz Analysis

Number of images: 8

Maximum 2D CAX->BB distance: 1.20mm

Median 2D CAX->BB distance: 0.86mm

Shift to iso: facing gantry, move BB: RIGHT 0.23mm; IN 0.30mm; DOWN 0.53mm

Gantry 3D isocenter diameter: 1.21mm (4/8 images considered)

Maximum Gantry RMS deviation (mm): 0.80mm

Maximum EPID RMS deviation (mm): 2.69mm








Gantry+Collimator 3D isocenter diameter: 1.49mm (6/8 images considered)

Collimator 2D isocenter diameter: 0.28mm (3/8 images considered)

Maximum Collimator RMS deviation (mm): 1.18

Couch 2D isocenter diameter: 0.93mm (3/8 images considered)

Maximum Couch RMS deviation (mm): 1.20

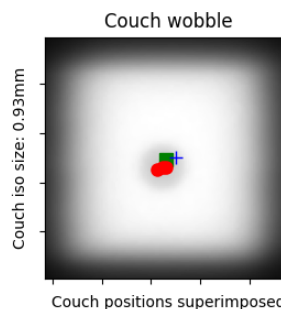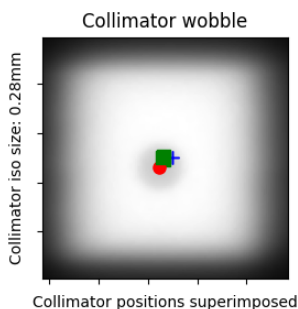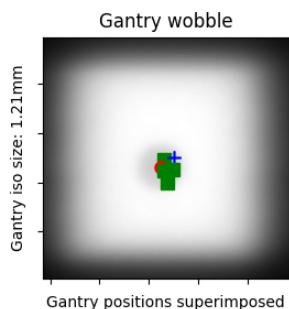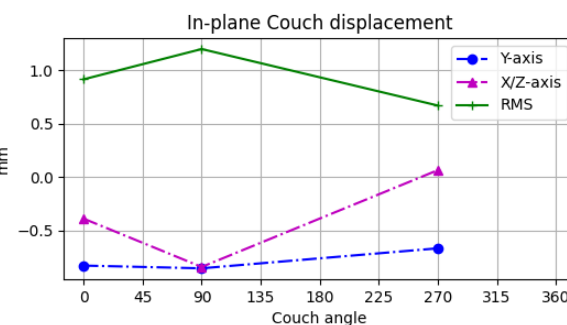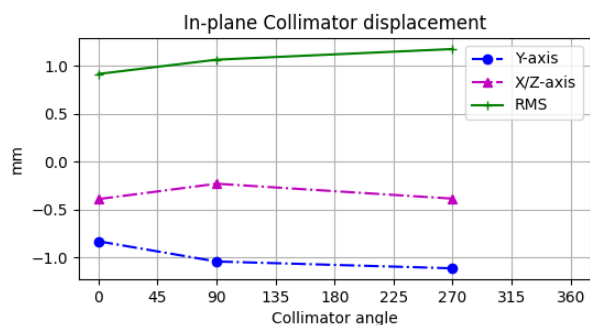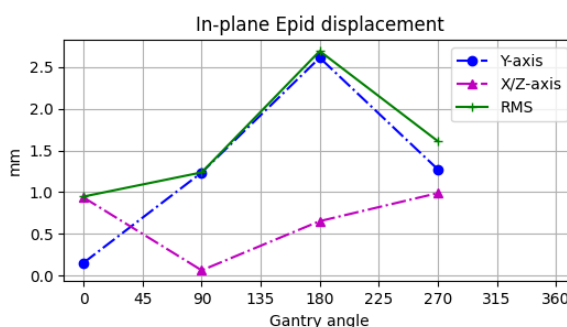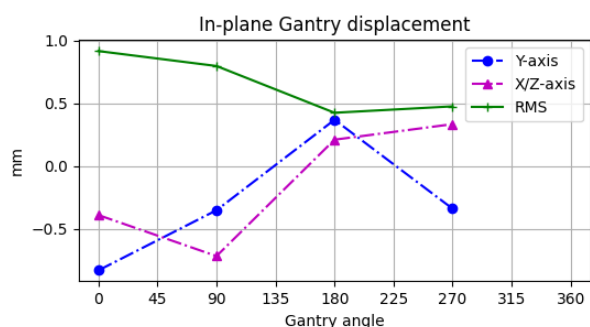

Notes:

Análise de Winston-Lutz com pylinac

Campo 2 cm X 2 cm

Tolerância: 1 mm

Limite de ação: 0.8 mm

Local das imagens: ./WL\WL_2019-11-28_analisado



Winston-Lutz Analysis

Metadata:


Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

Gantry images

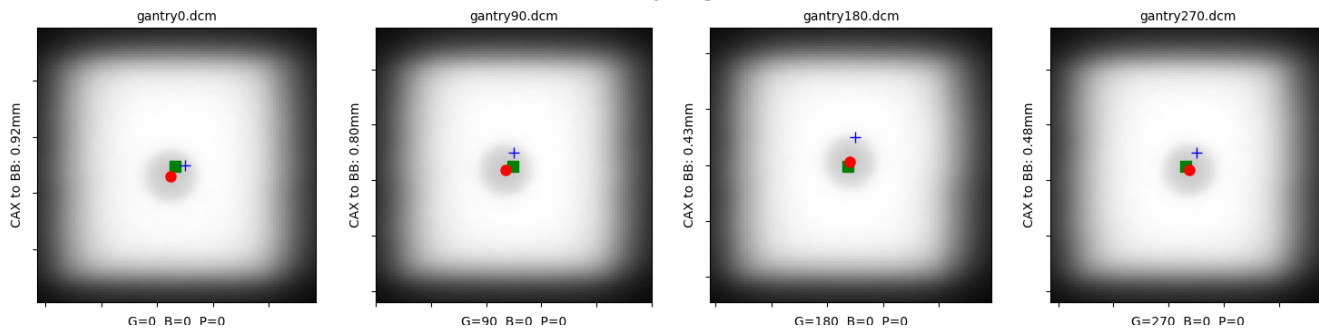



Winston-Lutz Analysis

Metadata:


Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

Collimator images

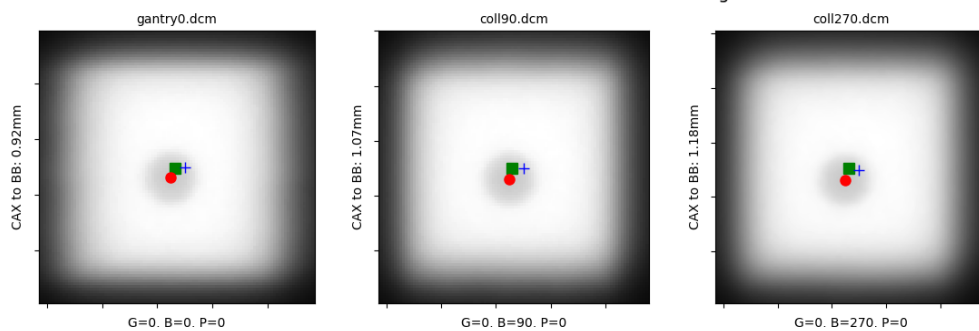



Winston-Lutz Analysis

Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

Couch images

