



Winston-Lutz Analysis

Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

Winston-Lutz Analysis

Number of images: 16

Maximum 2D CAX->BB distance: 1.10mm

Median 2D CAX->BB distance: 0.38mm

Shift to iso: facing gantry, move BB: RIGHT 0.03mm; OUT 0.16mm; UP 0.05mm

Gantry 3D isocenter diameter: 1.44mm (8/16 images considered)

Maximum Gantry RMS deviation (mm): 1.10mm

Maximum EPID RMS deviation (mm): 7.19mm








Gantry+Collimator 3D isocenter diameter: 1.44mm (14/16 images considered)

Collimator 2D isocenter diameter: 0.71mm (5/16 images considered)

Maximum Collimator RMS deviation (mm): 0.36

Couch 2D isocenter diameter: 1.10mm (4/16 images considered)

Maximum Couch RMS deviation (mm): 0.82

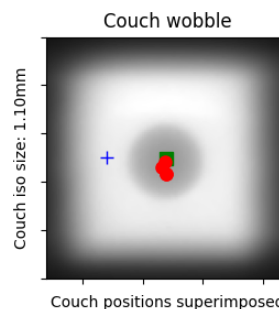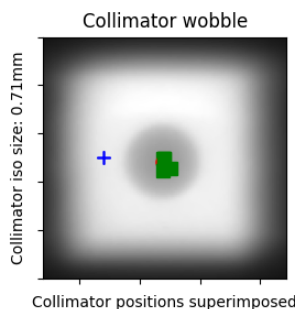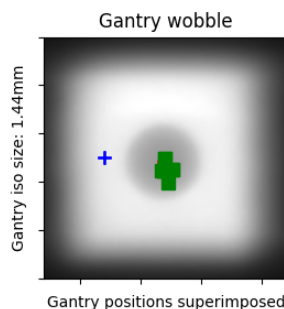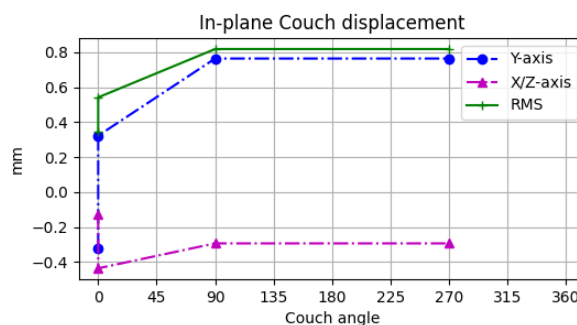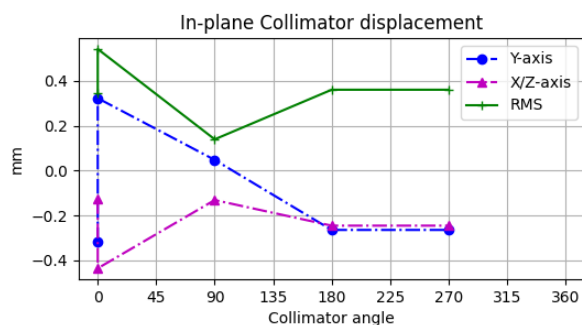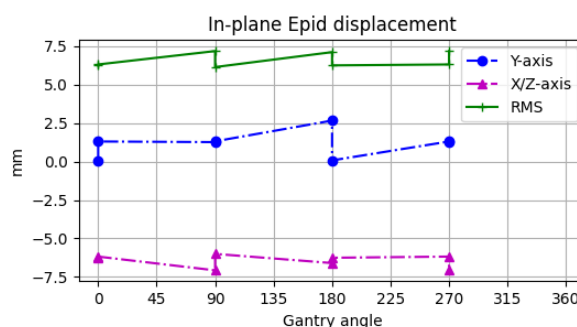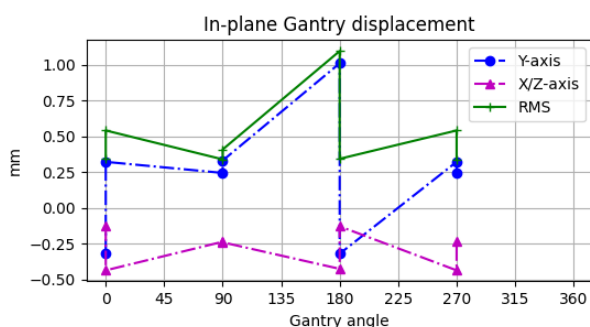

Notes:


Análise de Winston-Lutz com pylinac

Campo 2 cm X 2 cm

Tolerância: 1 mm

Limite de ação: 0.8 mm

Local das imagens: ./WL_IOSP_20201216\WL_2020-10-28_E

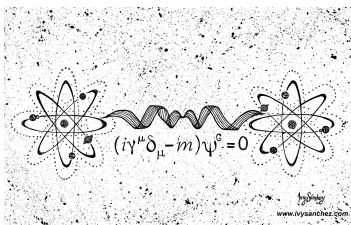

Winston-Lutz Analysis

Metadata:


Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

Gantry images

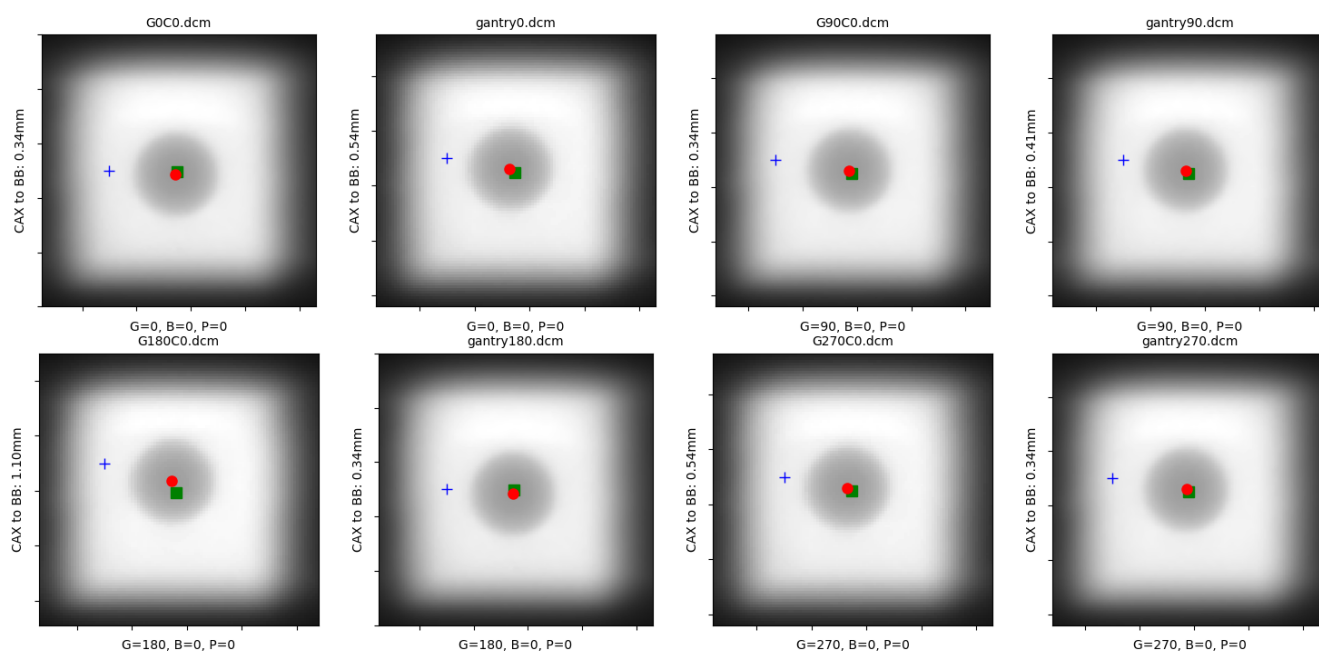



Winston-Lutz Analysis

Metadata:


Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

Collimator images

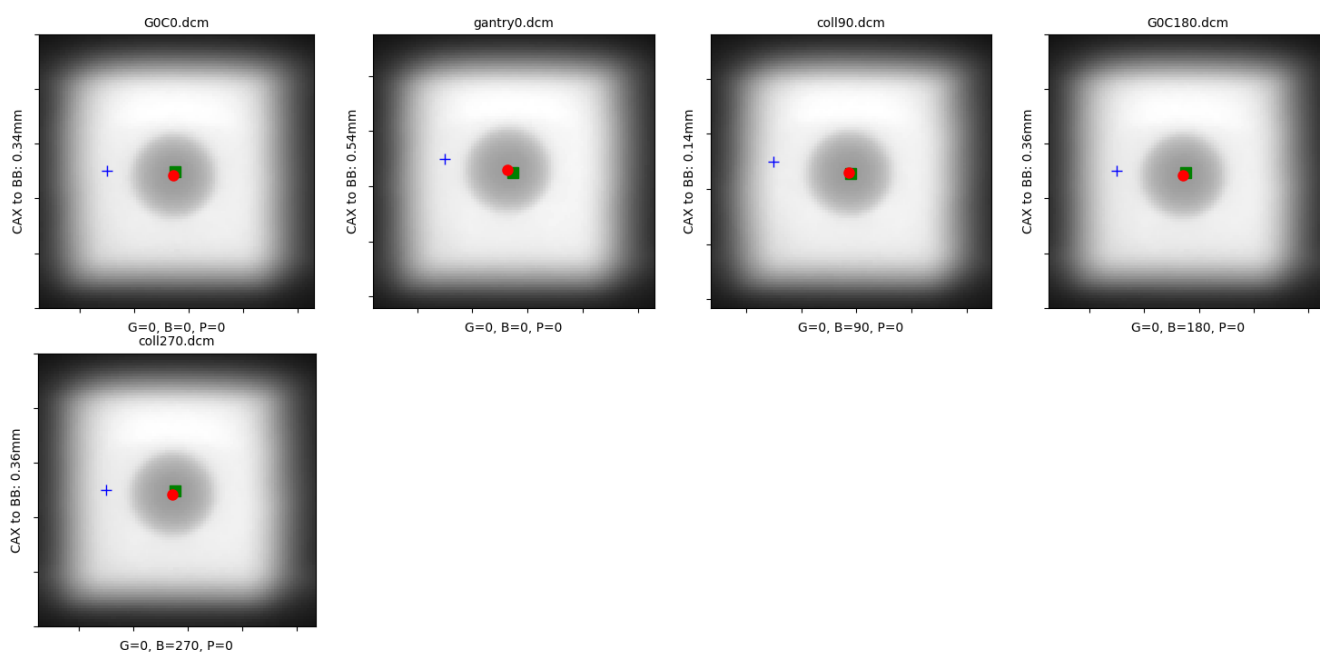



Winston-Lutz Analysis

Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

Couch images

